

PRIORITY		
_____	Mandatory	<u> x </u> As Required
_____	Next Visit	_____ At Installation
<u> x </u>	Information	

FIELD SERVICE BULLETIN

DATE ISSUED: September 4, 2014

NUMBER: 744

New Coherent GUI Versions for DPSS Glasgow and DPSS Santa Clara Products

PRODUCT DPSS Glasgow products and DPSS Santa Clara products

PURPOSE To inform Field Service, Representatives, and Coherent Subsidiaries that there are now two separate versions of the Coherent GUI.


DESCRIPTION A separate version of GUI for DPSS Glasgow products and DPSS Santa Clara products will allow both business units to make updates or alterations independently of each other without having to go through regression testing on each other's product variants.

Both applications can co-exist on the same PC. One is called **GLATransitionGUI** and the other is called **CoherentGUI** (Fig.1). Each version is available to download on two separate web pages for each application under the appropriate product family from the Coherent website (See following links).

[GLATransitionGUI](#) - (Chameleon/Talisker)

[Coherent GUI](#) - (Avia/Verdi/Vitesse/Datona)

Fig.1

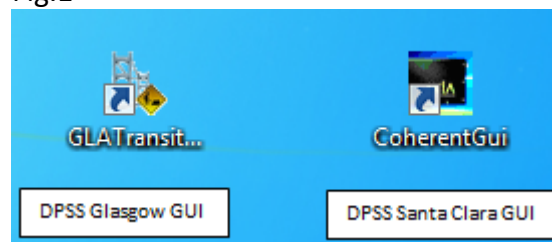

ADDITIONAL INFORMATION


Users of CoherentGUI v2.41 and below

1) Ignoring the update:

It is not essential to update if currently using the original CoherentGUI i.e. v2.41 and below (especially if you are interfacing DPSS Glasgow products), however you may be prompted to upgrade to the latest version each time the CoherentGUI is opened (notification shown in Fig.2).

In order to prevent this notification from reappearing it is advised to check the 'Don't check for updates when starting' box, then select 'Ignore this version'.

Fig.2 Ignoring the update

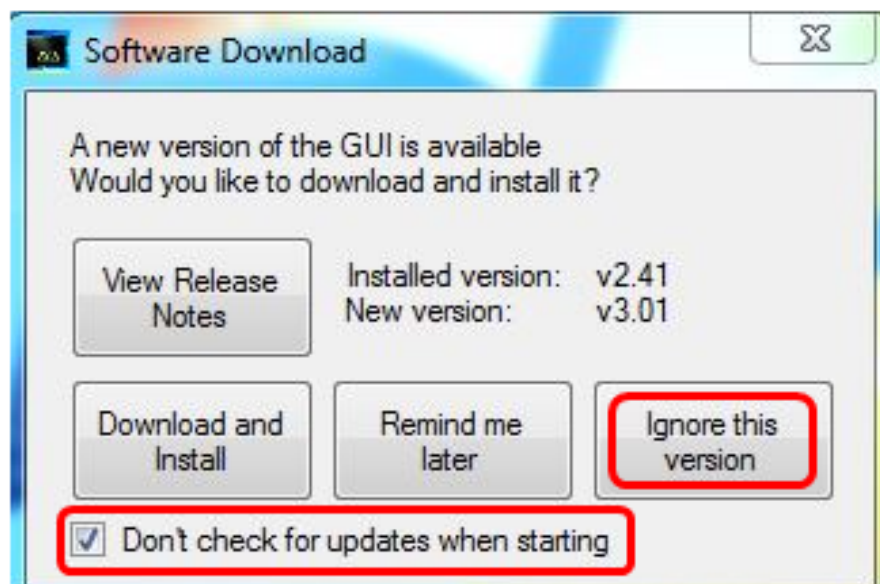

2) Installing the update:

By selecting 'Download and Install' (Fig.3) the CoherentGUI for DPSS Santa Clara products **only** will be installed. Once a Chameleon laser is connected to this GUI, then there will be an option to download the GLATransitionGUI for DPSS Glasgow products.

The following guide shows how to download both GUI versions by selecting the 'Download and Install' tab.

Fig.3 Select 'Download and Install'



Fig.4 CoherentGUI installer is downloading

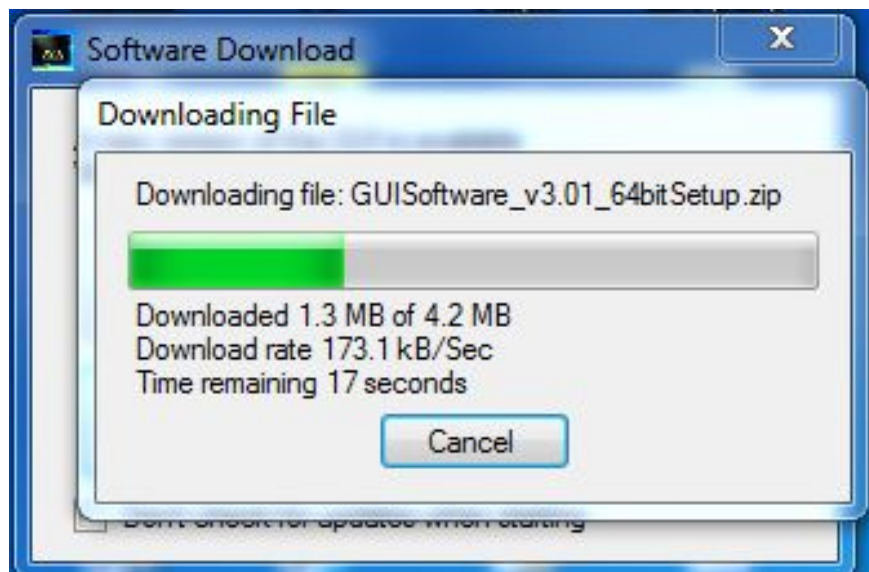


Fig.5 Select 'Next' to continue



Fig.6 Choose from two options then select 'Next'

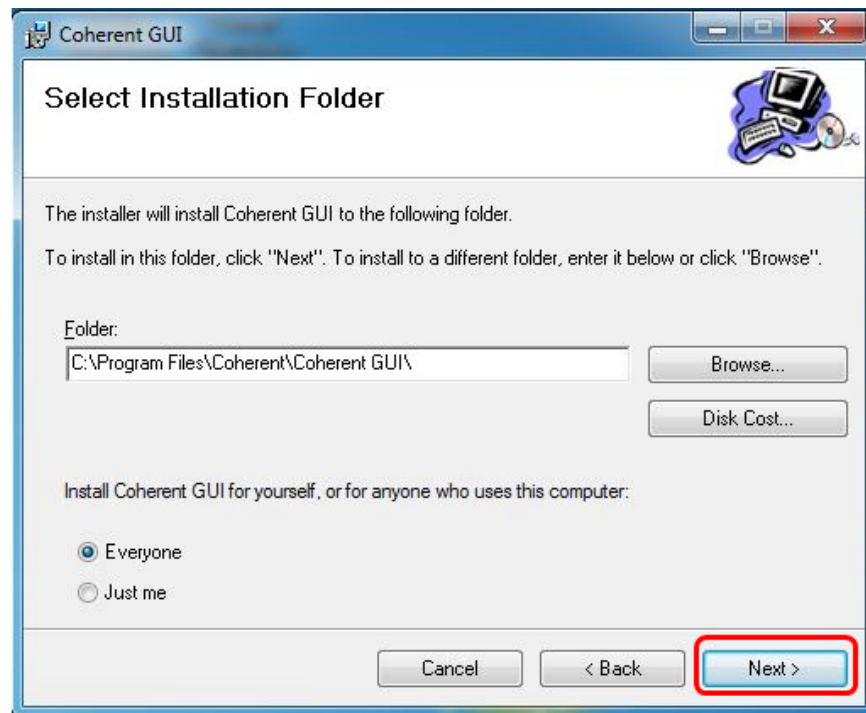


Fig.7 Select 'Next' to install the CoherentGUI on your computer

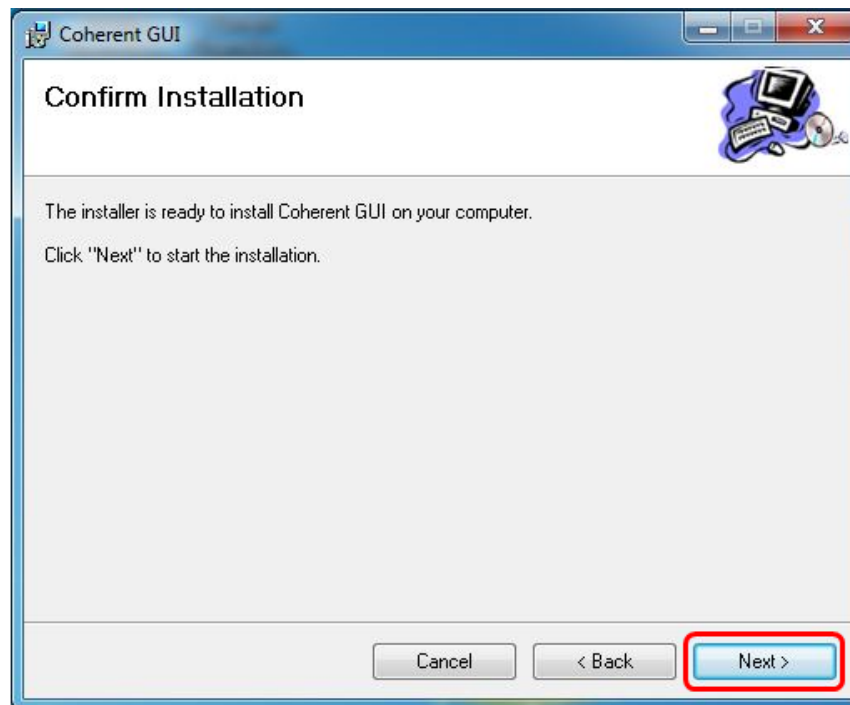


Fig.8 CoherentGUI is installing



Fig.9 The Installation is complete, select 'Close'

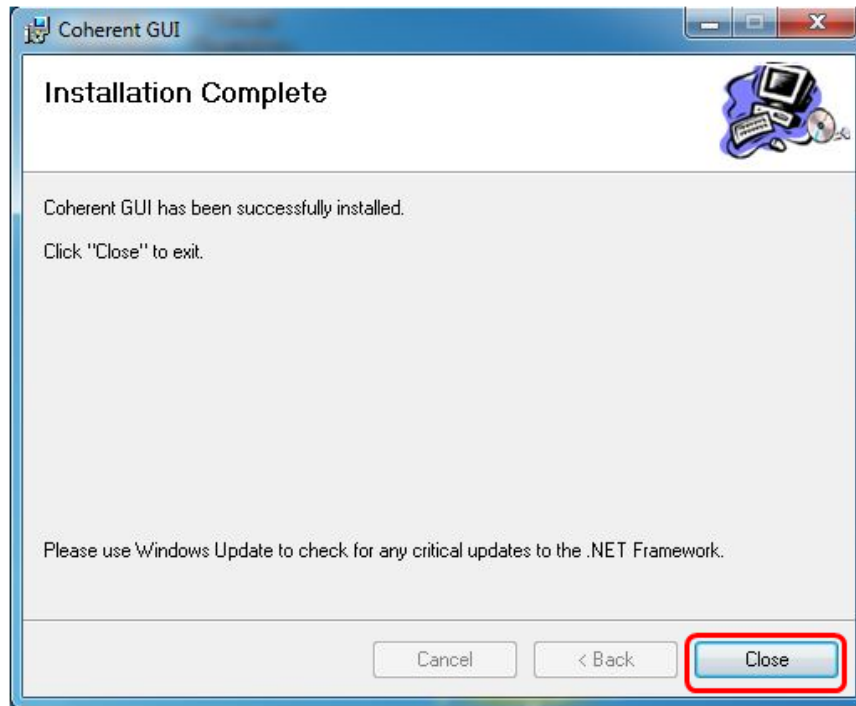


Fig.10 Double click the CoherentGUI icon on desktop



Fig.11 The CoherentGUI for DPSS Santa Clara products will open

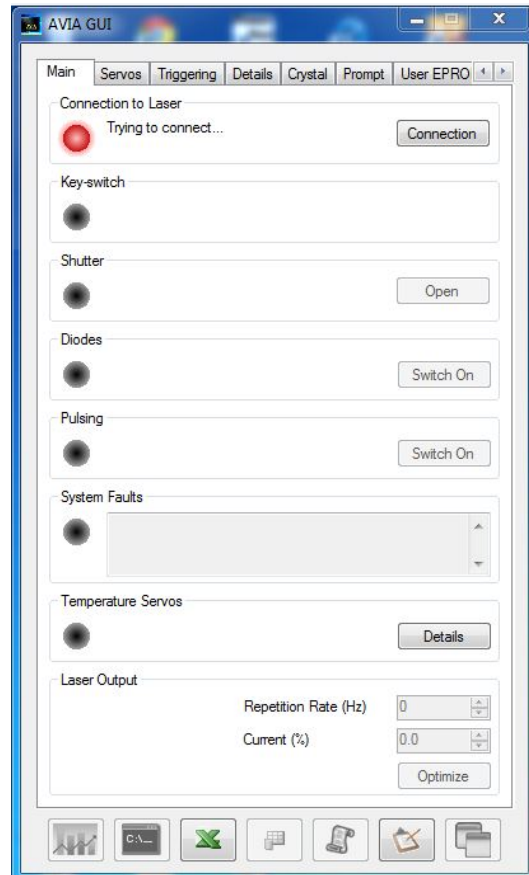


Fig.12 This notification will be displayed if a **Chameleon** laser is connected to the GUI. Select 'Yes' to go to the Chameleon website to download the GLATransitionGUI for DPSS Glasgow products

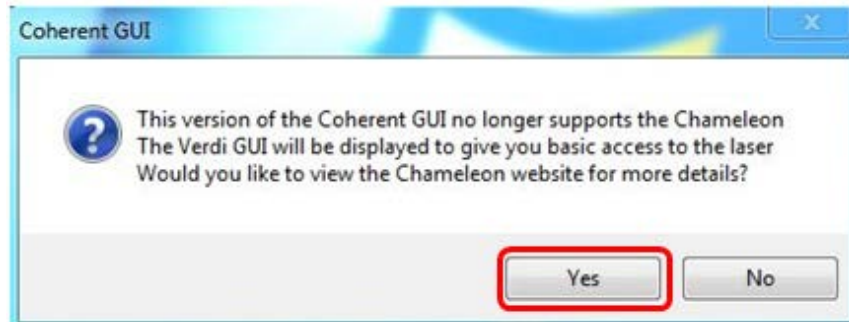


Fig.13 Select the GUI icon

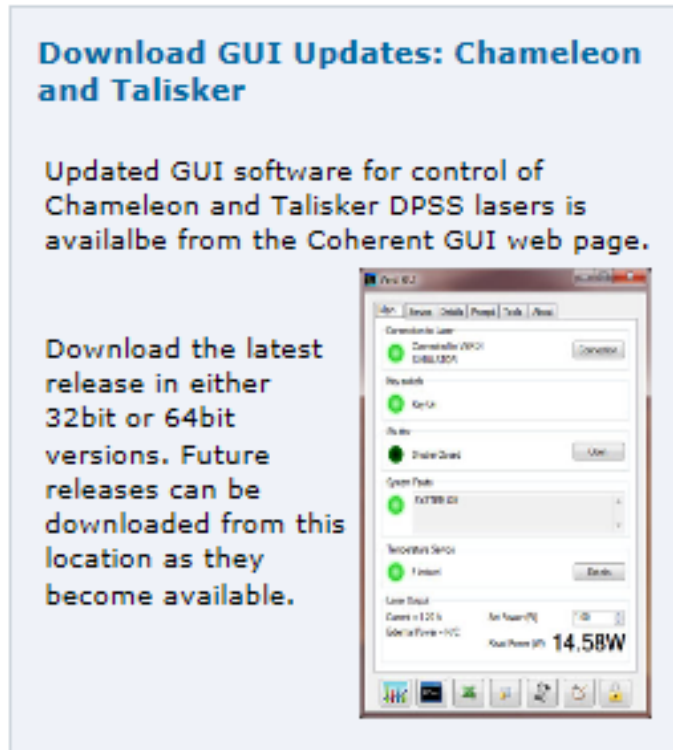


Fig.14 Choose the GUI S/W which is suitable for your PC to install the GLATransitionGUI for DPSS Glasgow products



PSE Initials

OJM

PRIORITY	
<input type="checkbox"/> Mandatory	<input type="checkbox"/> As Required
<input type="checkbox"/> Next Visit	<input type="checkbox"/> At Installation
<input checked="" type="checkbox"/> Information	

FIELD SERVICE BULLETIN

DATE ISSUED: April 25, 2014

NUMBER: 733


Coherent GUI v2.41 Error

DATE ISSUED 25th April 2014

PRODUCT Coherent GUI v2.41 32bit & 64bit

PURPOSE Inform Field Service Engineers of an error with the Coherent GUI v2.41 when used with Chameleon lasers systems. This occurs after the software is downloaded from the Coherent website.

DESCRIPTION 'The download of the JRE failed' appears while opening the Coherent GUI v2.41 after downloading the software from the Coherent website. The GUI can still be opened, however the spectrometer will not function and other random issues may occur.

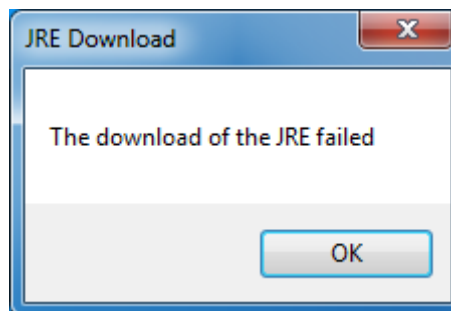

SOLUTION Use the 'GLA transition GUI' which can be found within Spectrum <https://spectrum.na.coherentinc.com/service/home/dpssScotlandSVC/SitePages/Home.aspx>.

For more information contact: PSE OM