

PRIORITY	
<input checked="" type="checkbox"/> Mandatory	<input type="checkbox"/> As Required
<input type="checkbox"/> Next Visit	<input checked="" type="checkbox"/> At Installation
<input checked="" type="checkbox"/> Information	

## FIELD SERVICE BULLETIN

**DATE ISSUED:** September 18, 2013

**NUMBER:** 710

### PolyScience Chillers

**PRODUCT** Evolution All

**PURPOSE** Information

**DESCRIPTION** As of September 2013, RLS will be switching from Lytron chillers to PolyScience chillers. Switching to these new chillers will require customers and FSEs to understand the differences between the Lytron and the PolyScience units. Most importantly, AC utility requirement changed to 208V/230V, 60Hz or 200V, 50Hz [115V units no longer applicable].

**ACTION** FSEs need to familiarize themselves with the PolyScience Operator's Manual, Manual Addendum and Service Procedures [see attached].

**PARTS** Table listed below.

**PSE Initials** KC

Old Chiller PN	Old Chiller Part Description	Old Chiller Where Used	Replacement Chiller PN	New Chiller Where Used	Note
1157141	Chiller, Air-Cooled Condenser, Single Phase, 825W 1.8gpm Flow @ 115VAC/60Hz, 690W 1.5gpm Flow @ 100VAC/50Hz	Evo-15 115V, 60hz or 100V, 50hz	1234960	All Evolutions 208V/230V, 60hz or 200V, 50hz	AC Utility Requirement Changed
1157142	Chiller, Air-Cooled Condenser, Single Phase, 690W 1.5gpm Flow @ 230VAC/50Hz	Evo-15 230V, 50hz	1234950	All Evolutions 240V, 50hz	
1157143	Chiller, Air-Cooled Condenser, Single Phase, 1050W 1.8gpm Flow @ 208/230VAC/60Hz, 875W 1.5gpm Flow @ 200VAC/50Hz	Evo-15 208V/230V, 60hz or 200V, 50hz	1234960	All Evolutions 208V/230V, 60hz or 200V, 50hz	
1157144	Chiller, Air-Cooled Condenser, Single Phase, 1650W 4.8gpm Flow @ 115VAC/60Hz, 1375W 4.0gpm Flow @ 100VAC/50Hz	Evo-HE 115V, 60hz or 100V, 50hz	1234960	All Evolutions 208V/230V, 60hz or 200V, 50hz	AC Utility Requirement Changed
1157145	Chiller, Air-Cooled Condenser, Single Phase, 1375W 4.0gpm Flow @ 230VAC/50Hz	Evo-HE 230V, 50hz	1234950	All Evolutions 240V, 50hz	
1157146	Chiller, Air-Cooled Condenser, Single Phase, 1650W 4.8gpm Flow @ 208/230VAC/60Hz, 1375W 4.0gpm Flow @ 200VAC/50Hz	Evo-HE 208V/230V, 60hz or 200V, 50hz	1234960	All Evolutions 208V/230V, 60hz or 200V, 50hz	
1157147	Chiller, Air-Cooled Condenser, Single Phase, 1650W 2.8gpm Flow @ 115VAC/60Hz, 1375W 2.3gpm Flow @ 100VAC/50Hz	Evo-30/45 115V, 60hz or 100V, 50hz	1234960	All Evolutions 208V/230V, 60hz or 200V, 50hz	AC Utility Requirement Changed
1157148	Chiller, Air-Cooled Condenser, Single Phase, 1375W 2.3gpm Flow @ 230VAC/50Hz	Evo-30/45 230V, 50hz	1234950	All Evolutions 240V, 50hz	
1157149	Chiller, Air-Cooled Condenser, Single Phase, 1650W 2.8gpm Flow @ 208/230VAC/60Hz, 1375W 2.3gpm Flow @ 200VAC/50Hz	Evo-45 208V/230V, 60hz or 200V, 50hz	1234960	All Evolutions 208V/230V, 60hz or 200V, 50hz	
1157150	Chiller, Water-Cooled Condenser, Single Phase, 1650W 4.8gpm Flow @ 115VAC/60Hz, 1375W 4.0gpm Flow @ 100VAC/50Hz	Evo-HE 115V, 60hz or 100V, 50hz	1235060	All Evolutions 208V/230V, 60hz or 200V, 50hz	AC Utility Requirement Changed
1157151	Chiller, Water-Cooled Condenser, Single Phase, 1375W 4.0gpm Flow @ 230VAC/50Hz	Evo-HE 230V, 50hz	1235050	All Evolutions 240V, 50hz	
1157152	Chiller, Water-Cooled Condenser, Single Phase, 1650W 4.8gpm Flow @ 208/230VAC/60Hz, 1375W 4.0gpm Flow @ 200VAC/50Hz	Evo-HE 208V/230V, 60hz or 200V, 50hz	1235060	All Evolutions 208V/230V, 60hz or 200V, 50hz	
1157153	Chiller, Water-Cooled Condenser, Single Phase, 1650W 2.8gpm Flow @ 115VAC/60Hz, 1375W 2.3gpm Flow @ 100VAC/50Hz	Evo-30/45 115V, 60hz or 100V, 50hz	1235060	208V/230V, 60hz or 200V, 50hz	AC Utility Requirement Changed
1157154	Chiller, Water-Cooled Condenser, Single Phase, 1375W 2.3gpm Flow @ 230VAC/50Hz	Evo-30/45 230V, 50hz	1235050	All Evolutions 240V, 50hz	
1157155	Chiller, Water-Cooled Condenser, Single Phase, 1650W 2.8gpm Flow @ 208/230VAC/60Hz, 1375W 2.3gpm Flow @ 200VAC/50Hz	Evo-45 208V/230V, 60hz or 200V, 50hz	1235060	All Evolutions 208V/230V, 60hz or 200V, 50hz	

---

# Operator's Manual

---

## Refrigerated Recirculating Chillers

---



# Table of Contents

<b>Introduction.....</b>	<b>3</b>
<b>General Information.....</b>	<b>4</b>
General Safety Information .....	4
Safety Recommendations .....	4
Unpacking Your Chiller .....	5
Regulatory and Compliance Testing.....	5
Contents.....	5
<b>Controls and Components .....</b>	<b>6</b>
<b>Quick Start.....</b>	<b>7</b>
<b>Installation and Startup.....</b>	<b>8</b>
Site Requirements.....	8
Ambient Temperature and Relative Humidity .....	8
Location .....	8
Clearance .....	8
Electrical Power .....	8
Optional Signal Inputs/Outputs .....	9
External Ambient Tracking Temperature Probe .....	9
RS232 / RS485 Serial Output .....	9
Remote I/O Port.....	9
Plumbing .....	9
Process Piping.....	9
Reservoir Drain.....	9
External Water Filter.....	9
Closed Sytem or Cooling Coil Setup.....	10
Open Bath System Setup .....	10
Startup.....	10
Process Coolant .....	10
Electrical Power.....	11
Starting Process Fluid Flow.....	11
<b>Normal Operation .....</b>	<b>13</b>
Selecting the Temperature Unit .....	13
Displaying and Adjusting the Set Point.....	13
Displaying and Adjusting the Ambient Tracking Offset.....	13
Displaying and Adjusting the Remote Control Temperature .....	14
Selecting the Pressure / Flow Rate Display and Units .....	14
Selecting the Internal / External Temperature Display .....	14
Setting Operational Parameters.....	15
High Temperature Limit (HL).....	16
Low Temperature Limit (LL) .....	16
High Ambient Temperature Limit (HA) .....	16
Maximum Fluid Pressure (FP).....	16
Minimum Flow Rate (FL) .....	17
Maximum External / Internal Temperature Differential (Sd).....	17
Auto-Refrigeration Temperature (AF).....	18
Specific Heat Capacity (SH).....	18
Remote Probe (rP) .....	18
Internal Calibration Offset (c1) .....	19
External Calibration Offset (c2) .....	19
Baud Rate (PC) .....	20
Diagnostic Probe Password (di).....	20
<b>Display, Alarm and Error Messages .....</b>	<b>20</b>
<b>Adjusting the High Pressure Bypass Setting .....</b>	<b>23</b>
<b>Calibrating the Flow Rate .....</b>	<b>23</b>
<b>Enabling / Disabling the Local Lockout .....</b>	<b>24</b>
<b>Routine Maintenance and Troubleshooting .....</b>	<b>25</b>
Routine Maintenance .....	25

Standard Magnetic Drive Centrifugal Pump .....	25
Turbine Pump .....	25
Positive Displacement Pump.....	25
Condenser, Air Vents and Reusable Filter .....	25
Fluid Filter.....	25
Fluid Level .....	25
Temperature Calibration.....	25
Troubleshooting .....	26
Diagnostic Mode .....	27
<b>Technical Information .....</b>	<b>28</b>
General Specifications (all Chillers).....	28
Pump Performance .....	28
Performance Specifications – 60Hz Chillers.....	29
1/4-HP, 1/3-HP and 1/2-HP Chillers.....	29
WhisperCool 3/4-HP and 1-HP Chillers .....	30
Performance Specifications – 50Hz Chillers.....	33
1/4-HP, 1/3-HP and 1/2-HP Chillers.....	33
WhisperCool 3/4 HP and 1 HP Chillers.....	35
<b>Diagrams and Schematics .....</b>	<b>37</b>
Electrical Wiring Diagram.....	37
Process Flow Schematic.....	38
<b>Replacement Parts .....</b>	<b>39</b>
<b>RS232 Communications .....</b>	<b>41</b>
<b>Equipment Disposal (WEEE Directive).....</b>	<b>42</b>
<b>Service and Technical Support .....</b>	<b>43</b>
<b>Warranty .....</b>	<b>43</b>

## Introduction

Your Recirculating Chiller provides cooling power for demanding applications and serves as an economical alternative to tap water cooling systems. Extremely easy to use and maintain, it combines technological innovation with precise temperature control to deliver reliable heat removal for a wide variety of applications.

Here are some of the features that make your Chiller so user-friendly:








- Microprocessor-based temperature controller
- Large, easy to read digital temperature display (°C or °F)
- One-touch temperature set point adjustment
- Digital pressure/flow rate display (PSI, kPa, GPM, LPM) with push-button selection
- Cool Command™ modulated refrigeration system for enhanced temperature stability and extended compressor life
- WhisperCool™ Environmental Control System (3/4 and 1 HP models) with variable speed fan to reduce operational noise and decrease energy consumption
- Centrifugal, positive displacement, or regenerative turbine pump

It will take you very little time to get your Recirculating Chiller installed and running. This manual is designed to guide you quickly through the process. We recommend that you read it thoroughly before you begin.

## General Information

### General Safety Information


When installed, operated, and maintained according to the directions in this manual and common safety procedures, your Chiller should provide safe and reliable heat removal. Please ensure that all individuals involved in the installation, operation, or maintenance of this unit read this manual thoroughly prior to working with the unit.

	This symbol alerts you to wide range of potential dangers.
	This symbol advises you of danger from electricity or electric shock.
	This symbol marks information that is particularly important.
	This symbol indicates alternating current.
	These symbols on the Power Switch / Circuit Breaker indicate that they place the main power supply ON / OFF.
	This symbol on the Power Switch indicates that it places the unit in a standby mode. It DOES NOT fully disconnect the unit from the power supply.
	This symbol indicates a protective conductor terminal.

**Read all instructions pertaining to safety, set-up, and operation.  
Proper operation and maintenance is the user's responsibility.**

### Safety Recommendations

To prevent injury to personnel and/or damage to property, always follow your workplace's safety procedures when operating this equipment. You should also comply with the following safety recommendations:

	<ul style="list-style-type: none"><li>• Always connect the power cord on this unit to a grounded (3-prong) power outlet. Make certain that the outlet is the same voltage and frequency as your unit.</li><li>• Never operate the unit with a damaged power cord.</li><li>• Always turn the unit OFF and disconnect Mains power before performing any maintenance or service.</li></ul>
---	---

## Unpacking Your Chiller

Your Chiller is shipped in a special carton. Retain the carton and all packing materials until the unit is completely assembled and working properly. Set up and run the unit immediately to confirm proper operation. Beyond one week, your unit may be warranty repaired, but not replaced. If the unit is damaged or does not operate properly, contact the transportation company, file a damage claim and contact the company where your unit was purchased immediately.



**CAUTION:** Keep unit upright when moving. Be sure to follow your company's procedures and practices regarding the safe lifting and relocation of heavy objects.

## Regulatory and Compliance Testing

### CSA UL (60Hz units)

CAN/CSA C22.2 No. 61010-1-04 — Safety Requirements for Electrical Equipment for Measurement, Control and Laboratory Use, Part I: General Requirements.

CAN/CSA C22.2 No. 61010-010-04 — Safety Requirements for Electrical Equipment for Measurement, Control and Laboratory Use - Part 2-010: Particular Requirements for Laboratory Equipment for the Heating of Materials.

UL Std No. 61010-1 — Electrical Equipment for Laboratory Use, Part I: General Requirements.

UL Std No. 61010A-2-010 — Electrical Equipment for Laboratory Use, Part 2: Particular Requirements for Laboratory Equipment for the Heating of Materials.

### CE (50Hz units)

EC Low Voltage Directive 2006/95/EC

EC Electromagnetic Compatibility Directive 2004/108/EC

IEC 61010-1-2001

IEC 61326:2005 / EN 61326 : 2006

### Highly Accelerated Life Test (HALT) and Vibration Tests per ASTM D4169-8 (All units)



## Contents

The following items have been included with your Chiller:

- Operator's Manual
- Warranty Card
- IEC Power Cord
- Two sets of Inlet/Outlet Adapters: ½ inch male NPT, 5/8 inch male NPT










# Controls and Components

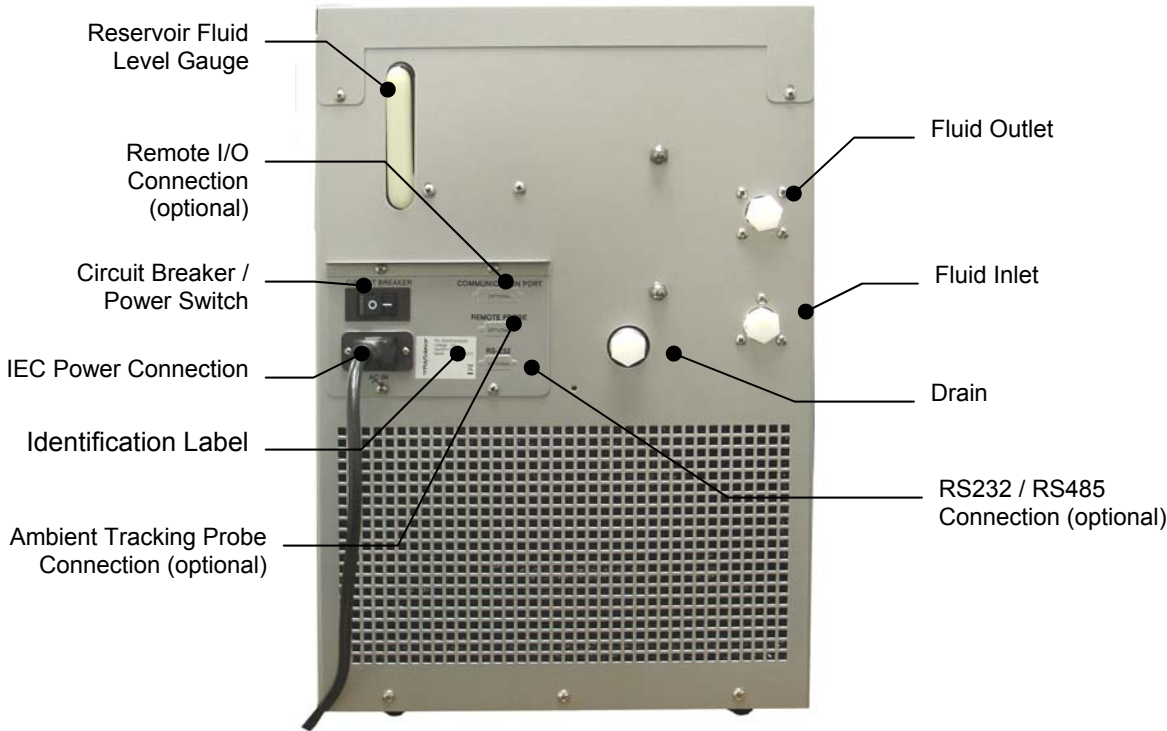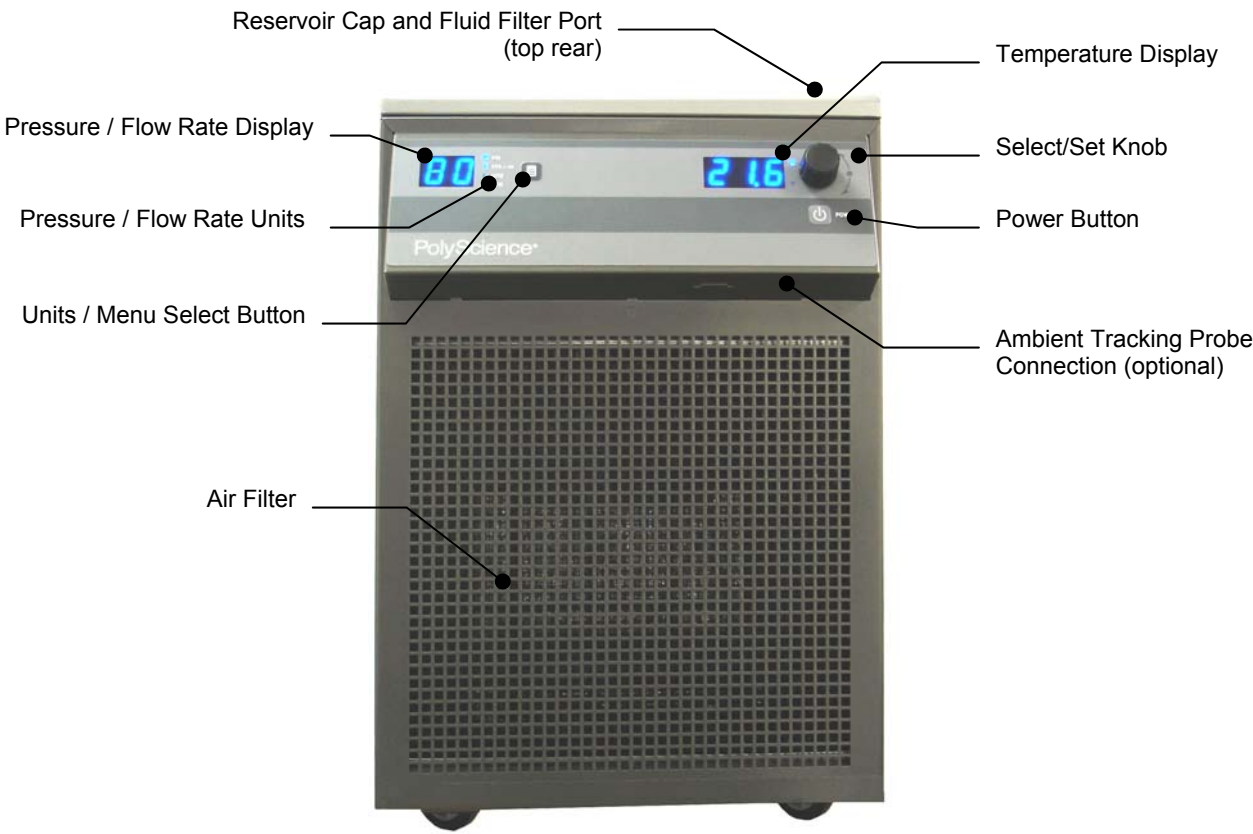

## Quick Start

See Installation & Startup for additional information.

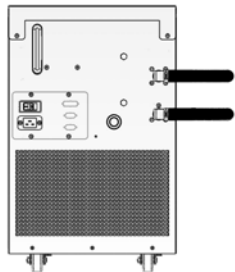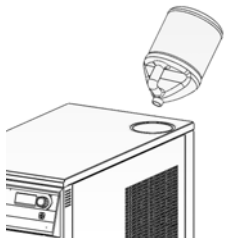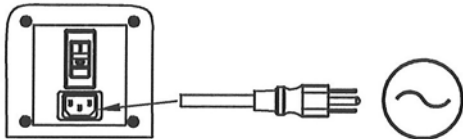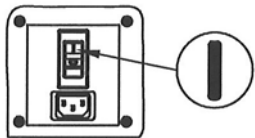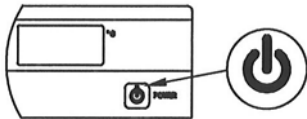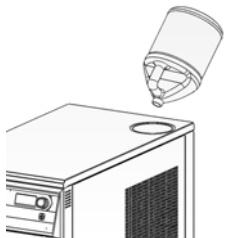1	Connect all process lines	
2	Fill reservoir with coolant	
3	Connect electrical power cord to Mains	
4	Turn Power Switch / Circuit Breaker ON	
5	Turn Controller ON	
6	Add coolant to reservoir as process lines fill	
7	Enter temperature set point	

## Installation and Startup



**WARNING:** Be sure all power is off before proceeding.

### Site Requirements

#### Ambient Temperature and Relative Humidity

The Chiller is designed for indoor installation in ambient temperatures between 5° and 30°C (41° and 86°F); relative humidity should not exceed 80% (non-condensing).

#### Location

The Chiller should be installed on a strong, level surface. It should be located as close to possible to the process requiring cooling. It should not be installed closer than 4 feet (1.4 meters) to a heat-generating source, such as heating pipes, boilers, etc. If possible, the Chiller should be located near a suitable drain to prevent flooding in the event of leaks. Do not place it where corrosive fumes, excessive moisture, excessive dust, or high room temperatures are present.

For ease of positioning and maneuverability, the Chiller is supplied with casters. The front wheels can be locked to keep the Chiller in place while in use.

To help prevent voltage drops, position the Chiller as close as possible to the power distribution panel. Avoid voltage drops by using a properly grounded power outlet wired with 14 gauge or larger diameter wire. The use of an extension cord is not recommended.



**NOTE:** The Chiller may be located at a level below that of the equipment being cooled. As long as the process remains closed, overflow will not occur when adding cooling fluid to the Chiller reservoir.

#### Clearance

Adequate clearance should be allowed on the front, sides, and rear of the Chiller for access to connections and components. The front and rear vents of the Chiller must be a minimum of 24 inches (61 cm) away from walls or vertical surfaces so air flow is not restricted.

#### Electrical Power

An IEC power cord is provided with the Chiller. It should be attached to the receptacle on the rear of the enclosure. Make sure that the power outlet used for the Chiller is properly grounded and matches the voltage and frequency indicated on the identification label on the back of the Chiller.

The use of an extension cord is not recommended. However, if one is necessary, it must be properly grounded and capable of handling the total wattage of the unit. The extension cord must not cause more than a 10% drop in voltage to the Chiller.



**WARNING:** DO NOT plug the Chiller into the electrical outlet until the unit is ready for startup (see *Startup*).

## Optional Signal Inputs/Outputs

### External Ambient Tracking Temperature Probe

This option allows you to control the cooling fluid temperature using an external temperature measurement (ambient room/machine temperature or process temperature). A 9-pin connector is provided on either the underside of the local control panel or on the rear panel for connecting the external probe.



**NOTE:** In order for the Chiller to properly recognize the presence of the external temperature probe, the probe must be connected to the unit before power is applied.

### RS232 / RS-485 Serial Output

This option allows you to remotely control the Chiller and/or output temperature readings to an external recorder or other auxiliary device. The maximum communications distance for Chillers equipped with the RS232 option is 50 feet (15 meters). The maximum distance for units equipped with the RS485 option is 4000 feet (1200 meters). A 9-pin D-connector is provided on the rear of the instrument enclosure for this connection.

### Remote I/O Port

This option allows you to use an external 12 VDC signal or a dry contact closure to turn the Chiller on and off. Chiller status is also available from this port. A 15-pin D-connector is provided on the rear of the instrument enclosure for this connection. See schematic at the end of this manual.

## Plumbing

### Process Piping

The Chiller has two internally threaded (1/2 inch ID NPT) fittings on the rear of the instrument housing for the process water connections. Two sets of adapters (1/2 inch ID and 5/8 inch ID) are supplied with the unit for connecting these fittings to the process piping.

To maintain a safe workplace and avoid leaks, special care should be taken when choosing hoses and connectors for the Chiller. It is the user's responsibility to ensure that the tubing and fittings connected to the Chiller are compatible with the fluid, temperature, and pressure being used.

**Pressure Ratings** — Hoses should be able to withstand the largest pressure that they will encounter. For "P" Series (positive displacement pump) and "T" Series (turbine pump) Chillers, this is 100 psi (689 kPa).

**Flexible Tubing** — Avoid tubing that will expand and take up fluid volume when operating at the desired pressure.

**Hose Diameter** — Process piping/hosing with a diameter smaller than 1/2 inch ID can be used if desired. However, keep in mind that using smaller diameter hosing increases pressure in the circulating system.

**Couplings and Clamps** — The use of screw-tightened hose clamps is necessary on all joints to insure good, tight connections. Quick connectors are not recommended as they have the potential for restricting flow rate.

### Reservoir Drain

A 1/2 inch NPT connection is provided for the reservoir's gravity drain. It should be piped to a drain or receptacle positioned below the bottom of the reservoir. If a receptacle is used, be sure it is of sufficient volume to hold all the water in the reservoir, process, and process lines.

### External Water Filter

An optional water filter is available that can be connected to the Chiller's fluid inlet or fluid outlet. Consult supplier for additional information.

## Closed System or Cooling Coil Setup

Connect the Chiller's inlet and outlet to the external apparatus with hoses or pipes. The direction of the flow through the system can be controlled by the way the connections are made. Fluid is drawn into the Chiller through the "Inlet" connection; fluid is pumped out of the Chiller through the "Outlet" connection.



**NOTE:** When Chillers with the standard magnetic drive centrifugal pump are connected to an external apparatus with a built-in shutoff, an external bypass loop assembly (Cat. No. 510-147) may be needed if operating below 20°C (68°F). This bypass assembly continues flow circulation to and from the pump even when the main flow to the external apparatus is blocked.

## Open Bath System Setup

Position the external tank at least two feet (0.6 meter) above the Chiller's inlet.

Install a shutoff valve on both the inlet and outlet of the Chiller. Place the valves in the closed position.

Connect the shutoff valves to the external tank using the tubing of equal diameter (1/2 inch minimum) and length. Use the same size fittings on both the inlet and outlet; this will ensure a balanced flow.

Cut the external end of the suction (inlet) tube into a "V" shape so that the tube will not seal itself against the wall of the external tank. Both the pressure and suction tubing should be securely fastened to the external tank to prevent movement during use. When using flexible tubing, the suction (inlet) tubing must have a wall thickness that will not collapse under vacuum, particularly when going around bends.

Fill the external bath (see *Startup, Process Coolant* for suitable fluids).

Fill the Chiller reservoir to the top of the filler neck and install the filter screen and cap. Tighten the cap until it is securely sealed.

## Startup

### Process Coolant

#### Suitable Fluids



**IMPORTANT:** Only use fluids that will satisfy safety, health, and equipment compatibility requirements. Caustic, corrosive, or flammable fluids must never be used.

The Chiller is designed to accommodate a variety of coolant fluids (water, glycol mixtures, etc). For most applications above 20°C (68°F), distilled water is satisfactory. For operation at or below 20°C (68°F), the Chiller must be protected with an antifreeze solution. Ethylene glycol (laboratory grade) and distilled water in a 50/50 mixture is satisfactory from +20° to -15°C (68° to 5°F). Select a fluid that is compatible with the Chiller's wetted parts (brass, stainless steel, polyethylene, EPDM rubber, and nylon).



**WARNING:** Do not use caustic, corrosive, or flammable fluids.



**WARNING:** Operation below 20°C (68°F) requires antifreeze in the circulation fluid.



**WARNING:** FOR CHILLERS WITH MAGNETIC DRIVE PUMPS ONLY: A low temperature fluid, such as a mixture of 50% ethylene or propylene glycol / 50% distilled water, or equivalent, must be used under all operating conditions.



**CAUTION:** Do not use the following fluids:

- Automotive antifreeze with additives\*\*
- Hard tap water\*\*
- Deionized water with a specific resistance > 1 meg ohm (except units with the DI water compatible plumbing)
- Any flammable fluids
- Concentrations of acids or bases
- Solutions with halides: chlorides, fluorides, bromides, iodides or sulfur
- Bleach (Sodium Hypochlorite)
- Solutions with chromates or chromium salts
- Glycerin
- Syltherm fluids
- \*\* At temperatures above 40°C, additives or mineral deposits can adhere to the heater. If deposits are allowed to build up, the heater may overheat and fail. Higher temperatures and higher concentrations of additives can hasten deposit build up.

### Filling the Reservoir

Remove the filler cap from the reservoir and, using a funnel, add fluid until it reaches the MAX line on the reservoir's fluid level gauge. Once the reservoir is full, remove the funnel but do not replace the cap at this time.

### Electrical Power

Plug the Chiller's power cord into an appropriate electrical outlet.

Place the Circuit Breaker/Power Switch on the rear of the instrument enclosure in the "On" position. Three decimal points will appear on the Temperature display; two decimal points will appear on the pressure/flow rate display.



### Starting Process Fluid Flow

Press the Power Button on the front panel. The system startup sequence will begin and proceed as follows:

1. The pump will turn on and fluid will begin circulating through the system. The set point temperature will appear briefly on the Temperature display; after a few seconds, it will be replaced by the actual fluid temperature. Fifteen to twenty seconds after power up, the compressor will begin operating.
2. Check for leaks.
3. With the pump running, the reservoir's fluid level will drop as the process and/or process cooling lines fill with fluid. Add fluid as follows:
4. **Closed Systems:** Slowly add fluid to the reservoir until the liquid level remains stable.
4. **Open Bath Systems:**
  - A. Set the flow alarm to at least 4 LPM (see "Minimum Flow Rate")
  - B. Open the inlet and outlet valves on the Chiller; the suction created by the pump should begin drawing fluid through the inlet tubing into the Chiller reservoir.
  - C. Once flow is established (no air bubbles in inlet tubing), close the inlet and outlet valves and turn the Chiller 'off'.
  - D. Remove the reservoir cap and check the level of the fluid in the reservoir. Add coolant until it is level with the top of the filler neck.



**CAUTION:** Always close the inlet and outlet valves before turning power to the Chiller 'off' or removing the reservoir cap to prevent the external reservoir from flooding the Chiller.

- E. Replace the reservoir cap, open the inlet and outlet valves, and restart the Chiller.
- F. Observe the liquid level in the external reservoir; adjust the valve on the Chiller outlet as required to maintain a stable fluid level.



**IMPORTANT:** When running an open loop system for extended periods, the fluid level in the Chiller reservoir should be checked periodically to avoid low fluid conditions.

To check the reservoir fluid level, close the inlet and outlet valves, turn the Chiller 'off', and remove the reservoir cap. Slowly open the inlet and outlet valves and allow fluid to drain from the external reservoir into the Chiller reservoir. Close the valves when the fluid level within the Chiller reservoir reaches the top of the filler neck. Add fluid to the external reservoir as required. Replace the reservoir cap, open the inlet and outlet valves, and turn the Chiller back on.

## Normal Operation



**SPECIAL FEATURE :** ¼ and 1 HP model Chillers are equipped with the WhisperCool™ Environmental Control System, which controls fan speed based on the heat load. You will notice the fan speed changing gradually during operation. This is especially beneficial in an environment where a low noise level is desirable.



**NOTE:** The Chiller incorporates a special “lockout” feature that can be enabled to prevent unauthorized or accidental set point and other operational changes. This feature is described in detail under “Enabling and Disabling the Local Lockout.” It should not be enabled until all operational parameters have been set.

### Selecting the Temperature Unit (°C or °F)

The LEDs adjacent to the Temperature Display indicate the unit (°C or °F) used for temperature displays. To change from °C to °F or vice versa, proceed as follows:

To change to °F — Place the Circuit Breaker/Power Switch on the rear of the instrument in the “Off” position. Press and hold the Units/Menu Select Button while returning the Circuit Breaker/Power Switch to the “On” position.

To change to °C — Place the Circuit Breaker/Power Switch on the rear of the instrument in the “Off” position. Press and hold the Power Button on the front panel while returning the Circuit Breaker/Power Switch to the “On” position.



**IMPORTANT:** All user settings, except baud rate and calibration offset, return to the original factory defaults when the temperature unit is changed. The Chiller’s temperature set point and various alarm settings should be reset to the desired values.

### Displaying and Adjusting the Set Point

Press the Select/Set Knob on the front panel. The current set point temperature will be displayed, and the decimal point at the bottom right of the display will flash, indicating that the temperature can be changed.

Rotate the Select/Set Knob until the desired set point temperature is displayed. The setting is accepted after the Set key is pressed a second time or automatically after a few seconds of inactivity.

When the set point is 20.0°C greater than the fluid temperature, the compressor and fan turn off. When the fluid temperature has risen to 5.0°C below the set point, the compressor and fan (enabled up to 50°C) will turn on and will remain on. If the set point is greater than the Auto-Refrigeration Temperature Set Point by 2°C, the compressor and fan will remain off. See “*Setting Operational Parameters, Auto-Refrigeration Temperature*”.



**NOTE:** Temperature set point cannot be displayed or changed when the optional ambient temperature tracking feature is installed and enabled. See “*Displaying and Adjusting the Ambient Tracking Offset*” below and “*Setting Operational Parameters, Remote Probe*”.

### Displaying and Adjusting the Ambient Tracking Offset



**NOTE:** Ambient tracking is an optional function that may or may not be available on your Chiller. It permits you to control fluid temperature based on room or machine temperature plus or minus a user-adjustable offset temperature.

When the optional ambient tracking probe is installed and enabled (AtC, see “*Setting Operational Parameters, Remote Probe*”), the ambient tracking offset rather than the set point temperature is displayed when the Select/Set Knob on the front panel is pressed.

To change the displayed offset value, rotate the Select/Set Knob until the desired offset value is displayed. An offset value from -5.0°C to +5.0°C (-9.0° to +9.0°F) may be entered. The setting is accepted after the Select/Set knob is pressed a second time or automatically after a few seconds of inactivity.



## Displaying and Adjusting the Remote Control Temperature



**NOTE:** Remote temperature control is an optional function that may or may not be available on your Chiller. It permits you to control cooling based on the temperature of an external process.

When the optional remote control external probe is installed and enabled (rP, see *Setting Operational Parameters, Remote Probe*), the external temperature set point is displayed when the Select/Set knob on the front panel is pressed.

To change the external temperature set point, press and then rotate the Select/Set knob until the desired set point temperature is displayed. The setting is accepted after the Select/Set knob is pressed a second time or automatically after a few seconds of inactivity.

## Selecting the Pressure / Flow Rate Display and Units

The Chiller can be set up to display either fluid pressure (in PSI or kPa) or if fitted with flow sensor flow rate in GPM (Option E) or lpm (Option M). Pressing the Units/Menu Select button briefly toggles through the available selections.



**NOTE:** If there is no flow sensor fitted, the display will show “- -” when Option E or Option M is selected.

After lpm (Option E) is displayed, the next press of the button will display P.1 on the left hand display and the internal temperature on the right hand display. Pressing again will display P.2 on the left hand display and remote temperature if the remote probe option has been ordered and is installed. If the remote probe option has been ordered, but the probe is not fitted then right hand display will show “- -”. If the remote probe option has not been ordered, the right hand display will show “n-A”.



**NOTE:** The flow rate readout is intended as a reference only. If accurate flow reading are required, an external flow meter is recommended.



**NOTE:** Metric pressure reading output is displayed in kPa and must be multiplied by 100 for Pa.

## Selecting the Internal / External Temperature Display



**NOTE:** This section applies only when the ambient tracking probe or remote temperature control is installed and enabled. It allows the user to check or continuously display either the Chiller's internal outlet fluid temperature or the external ambient/process temperature.

When the ambient tracking probe is selected (AtC, see “*Setting Operational Parameters, Remote Probe*”), the Chiller normally displays the internal outlet fluid temperature. To display the external ambient temperature, press and release the Units/Menu Select button until P2 appears on the pressure/flow rate display.

When the remote temperature control probe is selected (rPC, see “*Setting Operational Parameters, Remote Probe*”), the Chiller normally displays the external process temperature. To display the internal outlet fluid temperature, press and release the Units/Menu Select button until P1 appears on the pressure/flow rate display.



**NOTE:** P1 or P2 will remain on the pressure/flow rate display until the Units/Menu Select button is pressed and released. The displayed temperature will revert to the default condition (internal temperature for the ambient probe, external temperature for the remote temperature control probe).

## Setting Operational Parameters

The Chiller's various operational parameters, such as temperature, flow rate, and pressure alarm values, are all user-adjustable. They are accessed by pressing and holding the Units/Menu Button until HL appears on the pressure/flow rate display. Pressing and releasing the Units/Menu Button once HL appears allows you to scroll through the various parameters; rotating the Select/Set knob allows you to change the displayed setting. You can accept the displayed value by either pressing the Select/Set knob or allowing the display to timeout.

Menu Item	Description	Reset by Unit Change (°C -°F)	Choice/Ranges / Comments	Default Setting
HL	High Temperature Limit Alarm Set Point	Yes	+15°C to 71°C / 58° to 159°F	71°C 159°F
LL	Low Temperature Limit Alarm Set Point	Yes	-46°C to +14°C / -50° to 58°F	-46°C -50°F
HA	High Ambient Temperature Alarm Set Point	Yes	+30° to 71°C. Always displayed and set in °C.	50°C
FP w/psi LED lit	Maximum Fluid Pressure Alarm Set Point	Yes	40 to 100 PSI	100 PSI
FP w/kPa LED lit	Maximum Fluid Pressure Alarm Set Point	Yes	280 to 690 kPa	680 kPa
FL w/Option E LED lit	Minimum Flow Rate Alarm Set Point	Yes	0 or 0.8 to 2.0 GPM	0.0 GPM
FL w/Option M LED lit	Minimum Flow Rate Alarm Set Point	Yes	0 or 3 to 8 LPM	0.0 LPM
Sd	Establishes the cooling and heating rate by setting a desired differential between the measured internal and external temperature.	Yes	4°C to 20°C Always displayed and set in °C NOTE: Sd is displayed only if the optional remote temperature control probe is installed.	5°C
AF	Auto-Refrigeration Temperature Set Point	No	+20° to 50°C Always displayed/set in °C.	40°C
SH	Specific Heat Capacity	Yes	0.10 to 2.00 BTU / LB°F	0.80 BTU / LB°F
rP	External temperature probe status	No	n-A (the remote probe is not available) nEP (the remote probe is not plugged in) NO (the remote probe is disabled) AtC (ambient temperature control probe enabled) rPC (remote temperature control probe enabled)	n-A
°1	Internal Probe Calibration Offset	No	±2.9°C. Always displayed/set in °C. Special access procedure required. See <i>Internal Calibration Offset</i>	0.0°C
°2	External Probe Calibration Offset (only appears if a second probe is connected)	No	±2.9°C. Always displayed/set in °C. Special access procedure required. See <i>External Calibration Offset</i>	0.0°C
PC	RS232 Baud Rate	No	00, 24, 48, 96, 192. Represents baud rates of 0 (no communication), 2400, 4800, 9600, and 19200.	96
di	Diagnostic Mode	No	See <i>Diagnostic Mode</i>	00

You can adjust the following settings for your particular application or simply accept the default values.

**High Temperature Limit (HL)**

This menu item serves two functions. First, it establishes the maximum allowable set point temperature and thus helps prevent an operator from inadvertently selecting a temperature set point above a pre-established value. Secondly, it serves as a high temperature alarm, automatically activating both audio and visual alarm indicators if the measured fluid temperature reaches the HL setting. This also causes the compressor, heater, fan, and pump to turn off.

To change the high limit value, rotate the Select/Set knob until the desired value is displayed.

HL		71.0	
----	---	------	---

**Low Temperature Limit (LL)**

This menu item also serves a dual function. First, it establishes the minimum allowable set point temperature and thus helps prevent an operator from inadvertently selecting a temperature set point below a pre-established value. Secondly, it serves as a low temperature alarm, automatically activating both audio and visual alarm indicators if the measured fluid temperature drops to the LL setting. This also causes the compressor, heater, fan, and pump to turn off.

To change the low limit value, rotate the Select/Set knob until the desired value is displayed.



LL		-46.0	
----	---	-------	---

**High Ambient Temperature Limit (HA)**

	<b>NOTE:</b> This value is always set in °C.		
---	--	--	--

This menu item protects the Chiller from overheating due to a high ambient temperature. Should the ambient temperature rise above the limit value, the audio and visual alarms will activate, and the compressor, heater, fan, and pump will turn off.


To change the high ambient temperature value, rotate the Select/Set knob until the desired value is displayed on the temperature readout.

HA		50	
----	---	----	---

**Maximum Fluid Pressure (FP)**

This is the maximum allowable fluid pressure and can be set in either PSI or kPa (the LED adjacent to the display indicates the active unit of measure). Should the fluid pressure rise above the maximum fluid pressure value, the audio and visual alarms will activate and the compressor, heater, fan, and pump will turn off.

To change the fluid pressure limit value, rotate the Select/Set knob until the desired value is displayed on the temperature readout.

	<b>NOTE:</b> When FP first appears, the PSI LED will be lit. To view the FP value in kPa, press the Units/Menu Button again. The FP will remain on the display, and the kPa LED will light up.		
---	--	--	--

FP	PSI		100	
----	-----	---	-----	---

FP	kPa x 100		6.8	
----	-----------	---	-----	---



**NOTE:** Chillers with positive displacement and turbine pumps also incorporate a built-in safety that automatically maintains fluid pressure below a valve-regulated pressure value. It maintains this maximum outlet pressure by diverting the flow of process fluid to the reservoir (i.e., begin internally recirculating the fluid). A maximum pressure value is set at the factory, but is user-adjustable. See “*Adjusting the High Pressure Bypass Setting*” for information on changing the maximum outlet pressure value.

### Minimum Flow Rate (FL)



**NOTE:** If flow rates are below 1.5 GPM (5.7 LPM), this feature is not recommended as nuisance alarms may result.

This is the minimum allowable flow rate and can be set in either GPM or LPM (the LED adjacent to the display indicates the active unit of measure). Should the fluid flow rate drop below the minimum value, the audio and visual alarms will activate, and the compressor, heater, fan, and pump will turn off.

To change the minimum flow rate value, rotate the Select/Set Knob until the desired flow rate value is displayed on the temperature readout.

With FL set to “0”, the flow alarm is disabled, and the chiller will continue to operate with the output flow blocked.



**NOTE:** When FL first appears, the GPM LED will be lit. To view the FL value in LPM, press the Units/Menu Button again. The FL will remain on the display, and the LPM LED will light up.

FL	• Option E		0.0	
----	------------	---	-----	---

FL	• Option M		0.0	
----	------------	---	-----	---

### Maximum External / Internal Temperature Differential (Sd) (optional)




**NOTE:** This menu item (Sd) is only applicable when the remote temperature control probe is installed and enabled (see “Specific Heat Capacity”. Values are always displayed in °C.

This value helps establish the cooling/heating rate when the remote temperature control probe is being used. Maximum external/internal differential temperature values from 4° to 20°C may be entered. The higher the setting, the more rapidly the Chiller will achieve the external temperature set point. Low differential temperature settings minimize the amount of temperature overshoot/undershoot that occurs when the measured external temperature reaches the external set point temperature.

Sd			5.0	
----	--	---	-----	---

**Auto-Refrigeration Temperature (AF)**

	<b>NOTE:</b> This value is always displayed/set in °C.
---	--

This menu item allows you to select the temperature at which refrigeration is activated. When the set point exceeds the auto-refrigeration temperature by more than 1.0°C, the cooling and the fan are turned off.

To change the auto-refrigeration temperature, rotate the Select/Set knob until the desired value is displayed.


AF			40.0	
----	--	---	------	---

**Specific Heat Capacity (SH)**

This menu item allows you to select the specific heat capacity of the fluid in BTU/LB °F. This optimizes the control performance of the Chiller. To change the specific heat capacity, rotate the Select/Set knob until the desired value is displayed.



SH			0.8	
----	--	---	-----	---

**Remote Probe (rP)**

	<b>NOTE:</b> If an external temperature probe is not installed, n-A will be displayed when this menu item is selected.
--	--

This menu item allows you to enable/disable the Chiller’s optional ambient temperature probe (AtC) or remote temperature control probe (rPC).

If you wish to operate the Chiller using the ambient tracking probe, rotate the Select/Set knob until AtC is displayed. When this setting is selected, the effective set point for the cooling fluid temperature will be the temperature sensed by the ambient tracking probe (this may be room or machine temperature) plus or minus an offset specified by the user (see “*Displaying and Adjusting the Ambient Tracking Offset*”).

rP		AtC	
----	---	-----	---

If you wish to operate the Chiller using the remote temperature probe, rotate the Select/Set knob until rPC is displayed. When this setting is selected, the fluid temperature is controlled according to the temperature sensed by the external probe. The rate of cooling is controlled through the maximum differential temperature setting (Sd) (see “*Maximum External / Internal Temperature Differential*”).

rP		rPC	
----	---	-----	---

If you do not wish to operate the Chiller using either of these external probes, rotate the Select/Set knob until NO is displayed.

rP		NO	
----	---	----	---

### Internal Calibration Offset (°C1)



**IMPORTANT:** To prevent the operator from accidentally changing the calibration offset, a special keystroke sequence is required to access this function.

This menu item allows you to adjust the Chiller's internal temperature reading to match that of a traceable standard. It allows you to offset the displayed temperature value by as much as  $\pm 2.9^{\circ}\text{C}$ .



**NOTE:** Calibration offset values are always set and displayed in  $^{\circ}\text{C}$ .

1. Press and hold the Units/Menu Button until HL appears on the display.
2. Press and release the Units/Menu Button until rP appears on the display.
3. Press and hold the Units/Menu Button.
4. While holding the Units/Menu Button, press and release the Select/Set knob.
5. When CL1 appears on the temperature readout, release the Units/Menu Button. The current calibration offset value will appear on the temperature readout and alternate with the fluid temperature reading (enabling you to simultaneously adjust the offset and see the effect on the temperature).
6. Rotate the Select/Set knob until the desired calibration offset is displayed. Press the Select/Set knob or simply allow the display to time out to accept the displayed value.

°C1		0.0	
-----	--	-----	--

### External Calibration Offset (°C2)



**IMPORTANT:** To prevent the operator from accidentally changing the calibration offset, a special keystroke sequence is required to access this function.

This menu item allows you to adjust the Chiller's external temperature reading to match that of a traceable standard. It allows you to offset the displayed temperature value by as much as  $\pm 2.9^{\circ}\text{C}$ . It appears only if the external temperature probe is installed.



**NOTE:** Calibration offset values are always set and displayed in  $^{\circ}\text{C}$ .

1. Press and hold the Units/Menu Button until HL appears on the display.
2. Press and release the Units/Menu Button until °C1 appears on the display.
3. Press and hold the Units/Menu Button.
4. While holding the Units/Menu Button, press and release the Select/Set knob.
5. When C2 appears on the temperature readout, release the Units/Menu Button. The current calibration offset value will appear on the temperature readout and alternate with the fluid temperature reading (enabling you to simultaneously adjust the offset and see the effect on the temperature).
6. Rotate the Select/Set knob until the desired calibration offset is displayed. Press the Select/Set knob or simply allow the display to time out to accept the displayed value.

°C2		0.0	
-----	--	-----	--

## Baud Rate (PC)

This menu item allows you to establish the baud rate for serial communication. Allowable settings are 0 (no serial communication), 24 (2400 baud), 48 (4800 baud), 96 (9600 baud), 192 (19200 baud).

To change the displayed setting, rotate the Select/Set knob until the desired baud rate is displayed.

Press the Select/Set Knob or allow the display to time out to accept the displayed value.

PC		96	
----	---	----	---

## Diagnostic Mode Password (di)

This menu item allows access to the diagnostic mode, which displays important operational information that can aid in troubleshooting. By default, it is set to 00 (see “*Diagnostic Mode*” for access instructions).

di		00	
----	---	----	---

## Display, Alarm and Error Messages

When certain conditions are detected, a message code flashes on the display and the local audio alarm sounds. Depending on the nature of the condition, power to various systems components, such as the compressor, heater, fan, and pump, is removed. When the condition is rectified, push the front panel power button or turn the circuit breaker off then on to clear the fault or error.

Message Code	Description	Action Required
<b>EAF</b>	Rear panel high ambient temperature (select models only)	Warning - The ambient temperature is higher than the set ambient limit. The display alternates between EHA and the fluid temperature; the unit continues normal operation. If this lasts for more than 5 seconds, Fault 17 occurs.  Lower ambient temperature.
<b>E C</b>	External remote control active, Chiller in standby (for units with remote control by 12 VDC option)	Normal — Unit idle until remotely activated.
<b>E-C</b>	External remote control active, Chiller in standby  Only appears when Chiller is equipped for remote control using a dry contact	Normal — Unit idle until remotely activated.
<b>EFL</b>	Low fluid level warning (for units with optional float switch)	Warning — Fluid flow switch is open, indicating low fluid in the reservoir. An alarm will sound once every 8 seconds for 5 occurrences. If the fluid level has not been raised after 40 seconds, Fault 5 occurs, and the unit turns off.
<b>EHA</b>	Front panel high ambient temperature warning.	Warning - The ambient temperature is higher than the set ambient limit. The display alternates between EHA and the fluid temperature; the unit continues normal operation. If the ambient temperature stays over 5°C above the high ambient limit for more than 5 seconds, Fault 16 occurs.  Lower ambient temperature or raise temperature limit.

Message Code	Description	Action Required
<b>EHL</b>	High temperature set point warning	Warning — The temperature set point is higher than the high temperature limit value. The display alternates between EHL and the fluid temperature; the unit continues normal operation. The high temperature limit alarm will be activated if the fluid temperature exceeds the HL value for 25 consecutive seconds.  Lower temperature set point or increase high limit value.
<b>ELL</b>	Low temperature set point warning	Warning — The temperature set point is lower than the low temperature limit value. The display alternates between ELL and the fluid temperature; the unit continues normal operation. The low temperature limit alarm will be activated if the fluid temperature exceeds the LL value for 25 consecutive seconds.  Increase temperature set point or decrease low limit value.
<b>LLO</b>	Local Lockout	Normal — Indicates that Local Lockout feature (see <i>Enabling / Disabling the Local Lockout</i> ) is enabled. Appears momentarily when Select/Set knob is pressed to view/change set point value.
<b>CAn</b>	Cancel Local Lockout	Normal — Indicates the Local Lockout feature (see <i>Enabling / Disabling the Local Lockout</i> ) has been disabled. Appears momentarily when Local Lockout status is changed from enabled (LLO) to disabled.
<b>LO-H2O</b>	No fluid flow and no fluid pressure	Warning — LO is displayed alternately with H2O and fluid temperature if the fluid pressure exceeds the high pressure limit or if the fluid flow is less than the low flow limit. If either error lasts for more than 10 seconds, Fault 8 or Fault 7 (respectively) occurs.

If a fault should occur, the left display will show Ft, and the right display will show one of the fault codes shown below.


Fault Code	Description	Action Required
<b>01</b>	Factory reserved	None.
<b>02</b>	Low limit temperature alarm	Alarm — Process fluid temperature is below the low temperature limit value for more than 25 seconds. Compressor, heater, and fan are turned off; pump remains on.  To clear the fault, turn the unit off then on using the panel power button, and decrease the LL value.
<b>03</b>	High limit temperature alarm	Alarm — Process fluid temperature is above the high temperature limit value for more than 25 seconds. Compressor, heater, and fan are turned off; pump remains on.  To clear the fault, turn the unit off then on using the panel power button, and increase the HL value.
<b>04</b>	Over-temperature protection alarm	Alarm — Process fluid temperature is above Chiller's factory set high temperature safety cutoff. Heater, compressor, and fan turned off; pump remains on.  Lower process temperature.
<b>05</b>	Low liquid level alarm (select models only)	Delayed Alarm — Activated when the liquid level in the reservoir falls below an acceptable level for over 40 seconds. Compressor, heater, fan, and pump are turned off.  Add fluid to reservoir.




Fault Code	Description	Action Required
<b>06</b>	High bath temperature alarm	Alarm — Fluid temperature has exceeded 82°C (180°F) for more than 5 seconds. Compressor, heater, fan, and pump are turned off. Lower fluid temperature.
<b>07</b>	Low flow alarm	Alarm — Flow rate has dropped below minimum flow rate setting for more than 10 seconds. Compressor, heater, fan, and pump are turned off. Note: Disabled during first 2 minutes of operation. Correct cause of low flow rate or decrease minimum flow rate setting.
<b>08</b>	High pressure alarm	Alarm — Fluid outlet pressure has exceeded high-pressure limit value for over 10 seconds. Compressor, heater, fan, and pump are turned off. Decrease outlet pressure by removing blockage or increase high-pressure limit value.
<b>09</b>	System fault	Fault — One or more settings are out of range. Compressor, heater, fan, and pump are turned off. Contact service representative for corrective action.
<b>10</b>	Electronic power component fault (Triac)	Fault — Heater triac has failed for more than 10 seconds. Compressor, heater, fan, and pump are turned off. Contact supplier.
<b>11</b>	Internal probe fault	Fault — Main temperature control probe has failed for more than 4 seconds. Compressor, heater, fan, and pump are turned off. Contact supplier.
<b>12</b>	External temperature probe fault (select models only)	Fault — External temperature control probe has failed for more than 4 seconds. Compressor, heater, fan, and pump are turned off. Replace ambient tracking probe or operate instrument using internal temperature probe. Contact supplier if fault persists.
<b>13</b>	Communications fault	Fault — Internal electronics failure. Compressor, heater, fan, and pump are turned off. Contact supplier.
<b>14</b>	ADC fault, internal probe	Fault — ADC for internal temperature probe is reading faulty values. Compressor, heater, fan, and pump are turned off. Contact supplier.
<b>15</b>	ADC fault, external probe	Fault — ADC for external temperature probe is reading faulty values. Compressor, heater, fan, and pump are turned off. Contact supplier.
<b>16</b>	Front panel high ambient temperature alarm	Alarm — Ambient temperature at front panel has exceeded the high ambient temperature limit exceeds by more than 5°C for more than 5 seconds. Compressor, heater, fan, and pump are turned off. Lower ambient temperature or increase high ambient temperature limit value. See <i>High Ambient Temperature Limit</i>
<b>17</b>	Rear panel high ambient temperature alarm (select models only)	Alarm — Ambient temperature at rear panel has exceeded the factory-set high ambient limit for more than 5 seconds. Compressor, heater, fan, and pump turned off. Lower ambient temperature.

# Adjusting the High Pressure Bypass Setting

Chiller with positive displacement or turbine pumps incorporate an automatic safety to maintain outlet pressure below a valve-regulated pressure. This valve is adjustable and is located inside the Chiller housing.



**CAUTION:** There are exposed fan blades inside the Chiller housing. Exercise extreme care when accessing or adjusting any interior components.



**WARNING:** Hazardous voltages are present.

To access the high-pressure bypass valve, remove the two screws at the upper left and right corners of the Chiller's rear panel, slide the top panel back about 2-3 inches, and lift off. The regulator valve is located in the left rear corner of the unit.

The high-pressure bypass is adjusted as follows:

1. Set the low flow rate alarm value to zero (see *Setting Operational Parameters, Minimum Flow Rate*). This will prevent the unit from activating the flow alarm while you are adjusting the maximum pressure setting.
2. Completely block the Chiller's outlet flow. This should cause the outlet pressure to rise.
3. Set the Pressure/Flow Rate display to read either PSI or kPa.
4. Rotate the handle on the pressure valve until the desired maximum pressure setting is shown on the Pressure/Flow Rate display.
5. Reset the flow alarm value to the previous setting.
6. Return the Pressure/Flow Rate display to the previous setting.
7. Replace the top panel of the Chiller, being sure to secure the bayonet-style prongs on the front of the panel in the openings at the front of the unit. Reinsert the two screws that secure the top panel to the rear panel of the unit.

## Calibrating the Flow Rate

There may be minor differences between the Chiller's displayed flow rate and the actual flow rate as determined by a certified flow rate measurement device. The Chiller's displayed flow rate reading can be adjusted to match an external device by changing the flow rate gain coefficient setting as follows:

1. Press and hold the Units/Menu Button until HL appears on the display.
2. Press and release the Units/Menu Button until PC appears on the display.
3. Press and hold the Units/Menu Button; Di will appear on the pressure/flow rate display.
4. While holding the Units/Menu Button, press the Select/Set Knob. FAC will flash briefly on the temperature display and then 00 will appear.
5. Keep holding the Units/Menu Button and rotate the Select/Set Knob until 17 appears on the pressure/flow rate display and then press the Select/Set Knob. You are now in the Diagnostic Mode.
6. Press the Units/Menu Button until CP appears on the pressure/flow rate display.
7. Press and hold the Units/Menu Button; Fc will appear on the pressure/flow rate display.
8. Press the Select/Set Knob. CLF will flash briefly on the temperature readout and then the current gain coefficient value will appear.
9. Rotate the Select/Set Knob until the desired flow rate is displayed on the flow readout. Press the Select/Set Knob or simply allow the display to time out to accept the displayed value.

FC		6.28	
----	---	------	---

## Enabling / Disabling the Local Lockout

This feature is used to prevent unauthorized or accidental changes to set point and other operational values. When enabled, the values for the following functions can be displayed, but not changed:

- Temperature unit
- Temperature set point
- Ambient tracking offset
- Pressure / flow rate units

To enable the local lockout, press and hold the Select/Set Knob until LLO is displayed (approximately 5 seconds). Once enabled, LLO will appear momentarily when the Select/Set Knob is pressed to display the set point.

To disable the local lockout, press and hold the Select/Set knob until CAn appears momentarily as local lockout status changes from enabled (LLO) to disabled (approximately 5 seconds).




**IMPORTANT:** The Local Lockout feature does not prevent set point changes entered via the RS232 interface.

# Routine Maintenance and Troubleshooting

## Routine Maintenance

The Chiller is designed to require a minimum of periodic maintenance.

### Standard Magnetic Drive Centrifugal Pump

	<b>WARNING:</b> Hazardous voltages are present. Turn all power to the Chiller OFF and unplug the power cord from the electrical outlet.
---	---

When operating continuously, this pump should be oiled every six (6) months with SAE 20 oil. The pump incorporates two oil ports for this purpose.

To access the pump:

1. Turn both power switches off and unplug the power cord.
2. Remove the top panel of the housing (held in place with two bolts at the upper left and right corners of the rear panel).
3. Remove the housing's side panels by lifting them out of the housing frame.

### Turbine Pump

No lubrication required.


### Positive Displacement Pump

No lubrication required.

### Condenser, Air Vents and Reusable Filter

To keep the system operating at optimum cooling capacity, the condenser, the air vents, and reusable filter should be kept free of dust and dirt. They should be checked on a regular basis and cleaned as required.

The reusable filter is easily accessed from either the left or right side of the unit. Use a mild detergent and water solution to wash off any accumulated dust and dirt. Rinse and dry thoroughly before reinstalling.

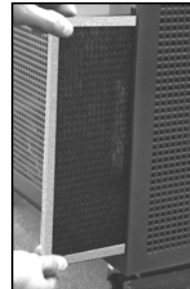

### Fluid Filter

A removable, highly efficient fluid filter is integrated into the fluid reservoir. To remove it for cleaning, simply remove the reservoir cap and lift the filter out of the reservoir. Rinse off accumulated particulate matter and reinstall.

### Fluid Level

The fluid level gauge on the rear of the Chiller should be periodically checked to determine if the fluid level needs to be topped off. Generally, fluid should be added whenever the level in the reservoir is at or near the "Low" gauge mark.

### Temperature Calibration

At times, there may be a minor temperature difference between the Chiller's displayed temperature and the actual temperature as determined by a certified temperature measurement device. There may also be situations where you want the displayed temperature to match a particular value to have standardization between different instruments. These adjustments can be performed using the Chiller's internal and/or external temperature calibration offset functions. See "*Setting Operational Parameters, Remote Probe and Internal Calibration Offset*".

## Troubleshooting



Many problems can be resolved by restoring the factory defaults. If this solves the problem, be careful when restoring your operational settings in order not to repeat the problem.

To restore the factory default settings:

1. Place the Power Switch/Circuit Breaker on the rear of the unit in the OFF position.
2. Press and hold the Power Button on the front panel while returning the Power Switch/Circuit Breaker to the ON position.



**WARNING:** Refer servicing to qualified service personnel.



**WARNING:** When electrical power is ON, dangerous voltages exist within chassis components. Use extreme care when measuring voltages on live circuits.

Problem	Possible Causes	Corrective Action
Unit does not run (digital displays blank)	No power to unit	Check that the electrical cord is secure and connected to an operating electrical outlet.  Check that Power Switch / Circuit Breaker on rear of unit is ON.
Unit does not run (three decimal points appear on temperature display, two decimal points on pressure/flow rate display)	Unit in Standby mode	Press Power Button on front panel.
No fluid circulation	Insufficient fluid in reservoir Blockage in circulating system Pump is not operating	Add fluid to reservoir. Remove blockage. Replace pump.
Insufficient circulation	Fluid viscosity too high External tubing diameter too small Restrictions in fluid lines Low line voltage	Replace with lower viscosity fluid. Replace with larger diameter tubing.  Check and correct as required. Check and correct as required.
Unit does not cool or cooling is insufficient	Dust build up on air filter or condenser Blocked air ventilation screens Excessive heat load  Ambient air temperature too high Low or high line voltage Faulty temperature sensor	Clean air filter and/or condenser as required.  Remove blockages as required. Check that heat load does not exceed capacity of chiller; correct as required. Decrease ambient air temperature.  Check and correct as required. Check the compressor upper, evaporator inlet and evaporator outlet temperature sensor readings (see "Diagnostic Mode"). If any of these temperature readings is -50°C, the sensor needs to be replaced.
Fault code 10 on display	Extreme electrical line interference  Triac failure	Plug unit into another power source. If problem persists, triac has failed.  Replace triac or triac driver as required.

Problem	Possible Causes	Corrective Action
Fault code 11 on display	Internal probe failure	Contact supplier
Fault code 12 on display	Loose external probe connection  Faulty external temperature probe	Check and correct as required.  Replace as required.  NOTE: Chiller may be operated using internal probe until problem is corrected.

## Diagnostic Mode

The Chiller incorporates a diagnostic mode, which displays important operational information that can aid in troubleshooting. To access the diagnostic mode, press and hold the Units/Menu button until HL appears on the Pressure/Flow Rate display. Scroll through the menu items by repeatedly pressing the Units/Menu button until the display reads PC. Press and hold the Units/Menu button until di appears on the display. Press and release knob. FAC should flash on the temperature readout.

The diagnostic menu now appears on the Pressure/Flow Rate display; the current value for the diagnostic item appears on the temperature readout. To toggle through the diagnostic menu items, press the Units/Menu button. To exit diagnostic mode, press the Select/Set knob.



**NOTE:** Diagnostic items are display values only; they cannot be changed..


Menu Item	Values	Description
Ut	xx	Compressor upper temperature (°C)
EI	xx	Evaporator inlet temperature (°C)
EO	xx	Evaporator outlet temperature (°C)
Li	xx	Secondary transformer voltage (V)
LF	xx	Line frequency (Hz)
SF	05	Fan control set point (°C)
EC	o00	External control (dry contact open or closed, voltage (0 - 36Vdc))
FF	ooo	Fluid status (fluid flow sensor, float switch(es))
At	xx	Ambient temperature reading (°C)
hP	36	Heat control "P" coefficient
cP	04	Cool control "P" coefficient
Fc	6.28	Flow gain
H%	xxx	Heater drive percentage
C%	xxx	PWM valve drive percentage
E%	xxx	Evaporator PID output
C%	xxx	Compressor PID output
F%	xxx	Fan drive percentage
Ct	xx	Chiller type (model)
Fb	hxx	Fuse bits (remote control voltage, contact closures, etc.)

# Technical Information


## General Specifications (all Chillers)

Temperature Set Point Resolution	±0.1°C
Temperature Stability	±0.1°C
Temperature Units	°C or °F
Pressure Units	PSI or kPa
Pressure Display Resolution Pressure Display Accuracy	1 PSI / 6.9 kPa ±3.5% of full scale (100PSI)
Flow Rate Units	GPM or LPM
Flow Rate Display Resolution Flow Rate Display Accuracy	0.1 GPM / 1 LPM +/- 0.4 GPM / 1.5 LPM
Pump Inlet and Outlet	½ inch NPT


## Pump Performance

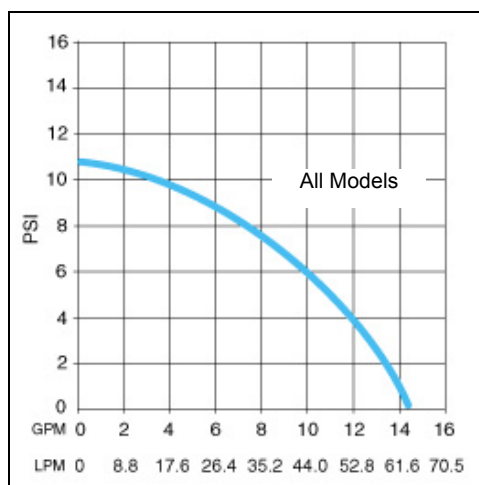

Magnetic Drive Centrifugal Pump

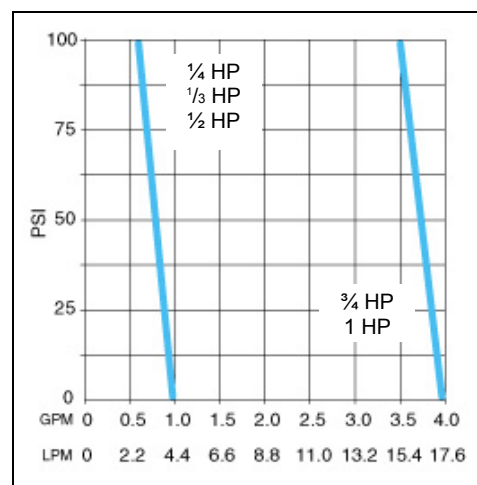

Positive Displacement Pump

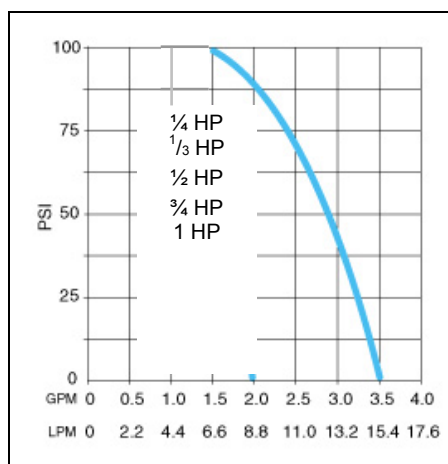

Turbine Pump

## Performance Specifications — 60Hz Chillers

### 1/4-HP, 1/3-HP and 1/2-HP Chillers

Magnetic Drive Centrifugal Pump (60Hz)						
Model:     Rfg = Refrigerating Only Rfg / Htg = Refrigerating & Heating						
Model Type	Rfg	Rfg / Htg	Rfg	Rfg / Htg	Rfg	Rfg / Htg
Operating Temperature	-10° to 40°C	-10° to 70°C	-10° to 40°C	-10° to 70°C	-10° to 40°C	-10° to 70°C
Cooling Temperature Range	-10° to 40°C	-10° to 50°C	-10° to 40°C	-10° to 50°C	-10° to 40°C	-10° to 50°C
Compressor	1/4 HP		1/3 HP		1/2 HP	
Cooling Capacity @ 20°C 10°C 0°C	950 watts 700 watts 300 watts	3244 BTU/hr 2391 BTU/hr 1025 BTU/hr	1430 watts 1000 watts 700 watts	4884 BTU/hr 3415 BTU/hr 2391 BTU/hr	1800 watts 1250 watts 850 watts	6147 BTU/hr 4269 BTU/hr 2903 BTU/hr
Pressure at 0 Flow Rate	10 psi / 69 kPa		10 psi / 69 kPa		10 psi / 69 kPa	
Flow Rate at 0 psi	4.1 gpm / 15.5 lpm		4.1 gpm / 15.5 lpm		4.1 gpm / 15.5 lpm	
Reservoir Capacity	1.1 gal / 4.2 liters		1.1 gal / 4.2 liters		1.1 gal / 4.2 liters	
Dimensions (h x w x d)	22-5/8 x 14-1/2 x 27-5/8 in. 57.5 x 36.8 x 70.2 cm.					
Shipping Weight	131 pounds 59.4 kg		143 pounds 64.8 kg		168 pounds 76.2 kg	
120 V, 60 Hz Chillers Volts Range	108 to 132V					
Amps	9.5A	10.0A	10.4A	10.7A	13.5A	13.8A



Positive Displacement Pump (60Hz)						
Model: Rfg = Refrigerating Only Rfg / Htg = Refrigerating & Heating						
Model Type	Rfg	Rfg / Htg	Rfg P	Rfg / Htg	Rfg	Rfg / Htg
Operating Temperature	-10° to 40°C	-10°C to 70°C	-10° to 40°C	-10°C to 70°C	-10° to 40°C	-10° to 70°C
Cooling Temperature Range	-10° to 40°C	-10°C to 50°C	-10° to 40°C	-10°C to 50°C	-10° to 40°C	-10° to 50°C
Compressor	1/4 HP		1/3 HP		1/2 HP	
Cooling Capacity @ 20°C 10°C 0°C	850 watts 600 watts 400 watts	2902 BTU/hr 2049 BTU/hr 1366 BTU/hr	1400 watts 990 watts 530 watts	4781 BTU/hr 3381 BTU/hr 1819 BTU/hr	1700 watts 1150 watts 750 watts	5806 BTU/hr 3927 BTU/hr 2561 BTU/hr
Flow Rate @ 0 psi	1 gpm / 3.8 lpm		1 gpm / 3.8 lpm		1 gpm / 3.8 lpm	
Pump Pressure (adjustable)	20 to 100 psi 138 to 689 kPa		20 to 100 psi 138 to 689 kPa		20 to 100 psi 138 to 689 kPa	
Reservoir Capacity	1.1 gal / 4.2 liters		1.1 gal / 4.2 liters		1.1 gal / 4.2 liters	
Shipping Weight	141 pounds 64 kg		153 pounds 69 kg		178 pounds 81 kg	
120 V, 60 Hz Chillers Volts Range	108 to 132V					
Amps	12.2A	12.5A	13.1A	13.5A	16.0A	16.4A

Turbine Pump (60Hz)						
Model: Rfg = Refrigerating Only Rfg / Htg = Refrigerating & Heating						
Model Type	Rfg	Rfg / Htg	Rfg	Rfg / Htg	Rfg	Rfg / Htg
Operating Temperature	-10° to 40°C	-10°C to 70°C	-10° to 40°C	-10°C to 70°C	-10° to 40°C	-10°C to 70°C
Cooling Temperature Range	-10° to 40°C	-10°C to 50°C	-10° to 40°C	-10°C to 50°C	-10° to 40°C	-10°C to 50°C
Compressor	1/4 HP		1/3 HP		1/2 HP	
Cooling Capacity @ 20°C 10°C 0°C	850 watts 600 watts 400 watts	2902 BTU/hr 2049 BTU/hr 1366 BTU/hr	1400 watts 990 watts 530 watts	4781 BTU/hr 3381 BTU/hr 1819 BTU/hr	1700 watts 1150 watts 750 watts	5806 BTU/hr 3927 BTU/hr 2561 BTU/hr
Flow Rate @ 0 psi	3.5 gpm / 13.2 lpm		3.5 gpm / 13.2 lpm		3.5 gpm / 13.2 lpm	
Pump Pressure (adjustable)	20 to 90 psi 138 to 621 kPa		20 to 90 psi 138 to 621 kPa		20 to 90 psi 138 to 621 kPa	
Reservoir Capacity	1.1 gal / 4.2 liters		1.1 gal / 4.2 liters		1.1 gal / 4.2 liters	
Shipping Weight	143 pounds 65 kg		156 pounds 71 kg		181 pounds 82 kg	
120 V, 60 Hz Chillers Volts Range	108 to 132V					
Amps	12.2A	12.5A	13.1A	13.5A	16.0A	16.4A

## WhisperCool 3/4-HP and 1-HP Chillers

Magnetic Drive Centrifugal Pump (60Hz)				
Model: Rfg = Refrigerating Only Rfg / Htg = Refrigerating & Heating				
Model Type	Rfg	Rfg / Htg	Rfg	Rfg / Htg
Operating Temperature	-10° to 40°C	-10° to 70°C	-10° to 40°C	-10 to 70°C
Compressor	3/4 HP		1 HP	
Cooling Capacity @ 20°C 10°C 0°C	2350 watts 1550 watts 975 watts	8026 BTU/hr 5294 BTU/hr 3330 BTU/hr	2900 watts 1835 watts 1100 watts	9904 BTU/hr 6267 BTU/hr 3757 BTU/hr
Pressure at 0 Flow Rate	10 psi / 69 kPa		10 psi / 69 kPa	
Flow Rate at 0 psi	4.1 gpm / 15.5 lpm		4.1 gpm / 15.5 lpm	
Reservoir Capacity	1.1 gal / 4.2 liters		1.1 gal / 4.2 liters	
Dimensions (h x w x d)	22-5/8 x 14-1/2 x 27-5/8 in. 57.6 x 36.8 x 70.2 cm.			
Shipping Weight	187 pounds 84.8 kg		189 pounds 85.7 kg	
208-230 V, 60 Hz Volts Range	187 to 253V			
Amps	9.2A	9.5A	9.5A	9.8A

Positive Displacement Pump (60Hz)				
Model: Rfg = Refrigerating Only Rfg / Htg = Refrigerating & Heating				
Model Type	Rfg	Rfg / Htg	Rfg	Rfg / Htg
Operating Temperature	-10° to 40°C	-10° to 70°C	-10° to 40°C	-10° to 70°C
Compressor	3/4 HP		1 HP	
Cooling Capacity @ 20°C 10°C 0°C	2300 watts 1550 watts 875 watts	7855 BTU/hr 5294 BTU/hr 2988 BTU/hr	2900 watts 1925 watts 1000 watts	9904 BTU/hr 6574 BTU/hr 3415 BTU/hr
Flow Rate @ 0 psi	3.5 gpm / 13.2 lpm		3.5 gpm / 13.2 lpm	
Pump Pressure (adjustable)	20 to 100 psi 138 to 689 kPa		20 to 100 psi 138 to 689 kPa	
Reservoir Capacity	1.1 gal / 4.2 liters		1.1 gal / 4.2 liters	
Shipping Weight	197 pounds 89 kg		199 pounds 90 kg	
208-230 V, 60 Hz Volts Range  Amps	187 to 253V			
	11.9A	12.2A	12.2A	12.5A

Turbine Pump (60Hz)				
Model: Rfg = Refrigerating Only Rfg / Htg = Refrigerating & Heating				
Model Type	Rfg	Rfg / Htg	Rfg	Rfg / Htg
Operating Temperature	-10° to 40°C	-10° to 70°C	-10° to 40°C	-10° to 70°C
Compressor	3/4 HP		1 HP	
Cooling Capacity @ 20°C 10°C 0°C	2300 watts 1550 watts 875 watts	7855 BTU/hr 5294 BTU/hr 2988 BTU/hr	2900 watts 1925 watts 1000 watts	9904 BTU/hr 6574 BTU/hr 3415 BTU/hr
Flow Rate @ 0 psi	3.5 gpm / 13.2 lpm		3.5 gpm / 13.2 lpm	
Pump Pressure (adjustable)	20 to 90 psi 138 to 621 kPa		20 to 90 psi 138 to 621 kPa	
Reservoir Capacity	1.1 gal / 4.2 liters		1.1 gal / 4.2 liters	
Shipping Weight	197 pounds 89 kg		199 pounds 90 kg	
208-230 V, 60 Hz Volts Range	187 to 253V			
Amps	11.9A	12.2A	12.2A	12.5A

Specifications subject to change without notice.

**Notes:** Refer to the serial number plate on the rear of the Chiller for model and electrical data.  
Cooling capacity (watts x 3.41) = BTU/hour. Performance specifications determined at ambient temperature of 20°C (68°F).  
Positive Displacement Pump Models: External pressure reducing assembly (Cat. No. 060302) steps down high outlet pressure to 10 to 45psi.

Environmental Conditions      Indoor use only  
Maximum Altitude:      2000 meters  
Operating Ambient:      5° to 30°C  
Relative Humidity:      80% for temperatures to 30°C  
Installation Category      II  
Pollution Degree:      2

## Performance Specifications — 50Hz Chillers

### 1/4-HP, 1/3-HP and 1/2-HP Chillers

Magnetic Drive Centrifugal Pump (50Hz)						
Model:     Rfg = Refrigerating Only Rfg / Htg = Refrigerating & Heating						
Model Type	Rfg	Rfg / Htg	Rfg	Rfg / Htg	Rfg	Rfg / Htg
Operating Temperature	-10° to 40°C	-10° to 70°C	-10° to 40°C	-10° to 70°C	-10° to 40°C	-10° to 70°C
Cooling Temperature Range	-10° to 40°C	-10° to 50°C	-10° to 40°C	-10° to 50°C	-10° to 40°C	-10° to 50°C
Compressor	1/4 HP		1/3 HP		1/2 HP	
Cooling Capacity @ 20°C 10°C 0°C	800 watts 600 watts 400 watts	2732 BTU/hr 2049 BTU/hr 1355 BTU/hr	1180 watts 960 watts 600 watts	4030 BTU/hr 3279 BTU/hr 2049 BTU/hr	1775 watts 1140 watts 765 watts	6062 BTU/hr 3893 BTU/hr 2613 BTU/hr
Pressure at 0 Flow Rate	9.5 psi / 66 kPa		9.5 psi / 66 kPa		9.5 psi / 66 kPa	
Flow Rate at 0 psi	3.9 gpm / 14.7 lpm		3.9 gpm / 14.7 lpm		3.9 gpm / 14.7 lpm	
Reservoir Capacity	1.1 gal / 4.2 liters		1.1 gal / 4.2 liters		1.1 gal / 4.2 liters	
Dimensions (h x w x d)	22-5/8 x 14-1/2 x 27-5/8 in. 57.5 x 36.8 x 70.2 cm.					
Shipping Weight	131 pounds 59.4 kg		143 pounds 64.8 kg		168 pounds 76.2 kg	
240 V, 50 Hz Chillers Volts Range Over Voltage	198 to 264V Category II					
Amps	5.6A	5.9A	5.9A	6.2A	7.2A	7.5A

Positive Displacement Pump (50Hz)						
Model: Rfg = Refrigerating Only Rfg / Htg = Refrigerating & Heating						
Model Type	Rfg	Rfg / Htg	Rfg P	Rfg / Htg	Rfg	Rfg / Htg
Operating Temperature	-10° to 40°C	-10°C to 70°C	-10° to 40°C	-10°C to 70°C	-10° to 40°C	-10° to 70°C
Cooling Temperature Range	-10° to 40°C	-10°C to 50°C	-10° to 40°C	-10°C to 50°C	-10° to 40°C	-10° to 50°C
Compressor	1/4 HP		1/3 HP		1/2 HP	
Cooling Capacity @ 20°C 10°C 0°C	700 watts 500 watts 300 watts	2391 BTU/hr 1708 BTU/hr 1025 BTU/hr	1280 watts 935 watts 485 watts	4371 BTU/hr 3193 BTU/hr 1656 BTU/hr	1775 watts 1140 watts 765 watts	6062 BTU/hr 3893 BTU/hr 2613 BTU/hr
Flow Rate @ 0 psi	0.95 gpm / 3.6 lpm		0.95 gpm / 3.6 lpm		0.95 gpm / 3.6 lpm	
Pump Pressure (adjustable)	20 to 95 psi 138 to 655 kPa		20 to 95 psi 138 to 655 kPa		20 to 95 psi 138 to 655 kPa	
Reservoir Capacity	1.1 gal / 4.2 liters		1.1 gal / 4.2 liters		1.1 gal / 4.2 liters	
Shipping Weight	141 pounds 64 kg		153 pounds 69 kg		178 pounds 81 kg	
240 V, 50 Hz Chillers Volts Range Over Voltage	198 to 264V Category II					
Amps						
	6.8A	7.1A	7.3A	7.6A	8.9A	9.2A

Turbine Pump (50Hz)						
Model: Rfg = Refrigerating Only Rfg / Htg = Refrigerating & Heating						
Model Type	Rfg	Rfg / Htg	Rfg	Rfg / Htg	Rfg	Rfg / Htg
Operating Temperature	-10° to 40°C	-10°C to 70°C	-10° to 40°C	-10°C to 70°C	-10° to 40°C	-10°C to 70°C
Cooling Temperature Range	-10° to 40°C	-10°C to 50°C	-10° to 40°C	-10°C to 50°C	-10° to 40°C	-10°C to 50°C
Compressor	1/4 HP		1/3 HP		1/2 HP	
Cooling Capacity @ 20°C 10°C 0°C	700 watts 500 watts 300 watts	2391 BTU/hr 1708 BTU/hr 1025 BTU/hr	1280 watts 935 watts 485 watts	4371 BTU/hr 3193 BTU/hr 1656 BTU/hr	1775 watts 1140 watts 765 watts	6062 BTU/hr 3893 BTU/hr 2613 BTU/hr
Flow Rate @ 0 psi	2.9 gpm / 11 lpm		2.9 gpm / 11 lpm		2.9 gpm / 11 lpm	
Pump Pressure (adjustable)	20 to 75 psi 138 to 572 kPa		20 to 75 psi 138 to 572 kPa		20 to 75 psi 138 to 572 kPa	
Reservoir Capacity	1.1 gal / 4.2 liters		1.1 gal / 4.2 liters		1.1 gal / 4.2 liters	
Shipping Weight	143 pounds 65 kg		156 pounds 71 kg		181 pounds 82 kg	
240 V, 50 Hz Chillers Volts Range Over Voltage	198 to 264V Category II					
Amps	6.8A	7.1A	7.3A	7.6A	8.9A	9.2A

## WhisperCool 3/4-HP and 1-HP Chillers

Magnetic Drive Centrifugal Pump (50Hz)				
Model: Rfg = Refrigerating Only Rfg / Htg = Refrigerating & Heating				
Model Type	Rfg	Rfg / Htg	Rfg	Rfg / Htg
Operating Temperature	-10° to 40°C	-10° to 70°C	-10° to 40°C	-10 to 70°C
Compressor	3/4 HP		1 HP	
Cooling Capacity @ 20°C 10°C 0°C	2250 watts 1600 watts 1075 watts	7684 BTU/hr 5464 BTU/hr 3671 BTU/hr	2750 watts 2050 watts 1400 watts	9238 BTU/hr 7001 BTU/hr 4781 BTU/hr
Pressure at 0 Flow Rate	9.5 psi / 66 kPa		9.5 psi / 66 kPa	
Flow Rate at 0 psi	3.9 gpm / 14.7 lpm		3.9 gpm / 14.7 lpm	
Reservoir Capacity	1.1 gal / 4.2 liters		1.1 gal / 4.2 liters	
Dimensions (h x w x d)	22-5/8 x 14-1/2 x 27-5/8 in. 57.6 x 36.8 x 70.2 cm.			
Shipping Weight	187 pounds 84.8 kg		189 pounds 85.7 kg	
240 V, 50 Hz Volts Range Over Voltage  Amps	198 to 264V Category II			
	9.2A	9.5A	9.5A	9.8A

Positive Displacement Pump (50Hz)				
Model: Rfg = Refrigerating Only Rfg / Htg = Refrigerating & Heating				
Model Type	Rfg	Rfg / Htg	Rfg	Rfg / Htg
Operating Temperature	-10° to 40°C	-10° to 70°C	-10° to 40°C	-10° to 70°C
Compressor	3/4 HP		1 HP	
Cooling Capacity @ 20°C 10°C 0°C	2200 watts 1500 watts 975 watts	7513 BTU/hr 5123 BTU/hr 3329 BTU/hr	2650 watts 1900 watts 1200 watts	9050 BTU/hr 6489 BTU/hr 4098 BTU/hr
Flow Rate @ 0 psi	2.9 gpm / 11 lpm		2.9 gpm / 11 lpm	
Pump Pressure (adjustable)	20 to 83 psi 138 to 572 kPa		20 to 83 psi 138 to 572 kPa	
Reservoir Capacity	1.1 gal / 4.2 liters		1.1 gal / 4.2 liters	
Shipping Weight	197 pounds 89 kg		199 pounds 90 kg	
240 V, 50 Hz Volts Range Over Voltage	198 to 264V Category II			
Amps	11.9A	12.2A	12.2A	12.5A

Turbine Pump (50Hz)				
Model: Rfg = Refrigerating Only Rfg / Htg = Refrigerating & Heating				
Model Type	Rfg	Rfg / Htg	Rfg	Rfg / Htg
Operating Temperature	-10° to 40°C	-10° to 70°C	-10° to 40°C	-10° to 70°C
Compressor	3/4 HP		1 HP	
Cooling Capacity @ 20°C 10°C 0°C	2200 watts 1500 watts 975 watts	7513 BTU/hr 5123 BTU/hr 3329 BTU/hr	2650 watts 1900 watts 1200 watts	9050 BTU/hr 6489 BTU/hr 4098 BTU/hr
Flow Rate @ 0 psi	2.9 gpm / 11 lpm		2.9 gpm / 11 lpm	
Pump Pressure (adjustable)	20 to 75 psi 138 to 572 kPa		20 to 75 psi 138 to 572 kPa	
Reservoir Capacity	1.1 gal / 4.2 liters		1.1 gal / 4.2 liters	
Shipping Weight	197 pounds 89 kg		199 pounds 90 kg	
240 V, 50 Hz Volts Range Over Voltage  Amps	198 to 264V Category II			
	11.9A	12.2A	12.2A	12.5A

Specifications subject to change without notice.

**Notes:** Refer to the serial number plate on the rear of the Chiller for model and electrical data.  
Cooling capacity (watts x 3.41) = BTU/hour. Performance specifications determined at ambient temperature of 20°C (68°F).  
Positive Displacement Pump Models: External pressure reducing assembly (Cat. No. 060302) steps down high outlet pressure to 10 to 45psi.


Environmental Conditions Indoor use only  
Maximum Altitude: 2000 meters  
Operating Ambient: 5° to 30°C  
Relative Humidity: 80% for temperatures to 30°C  
Installation Category: II  
Pollution Degree: 2

## Electrical Wiring Diagram





Process Flow Schematic

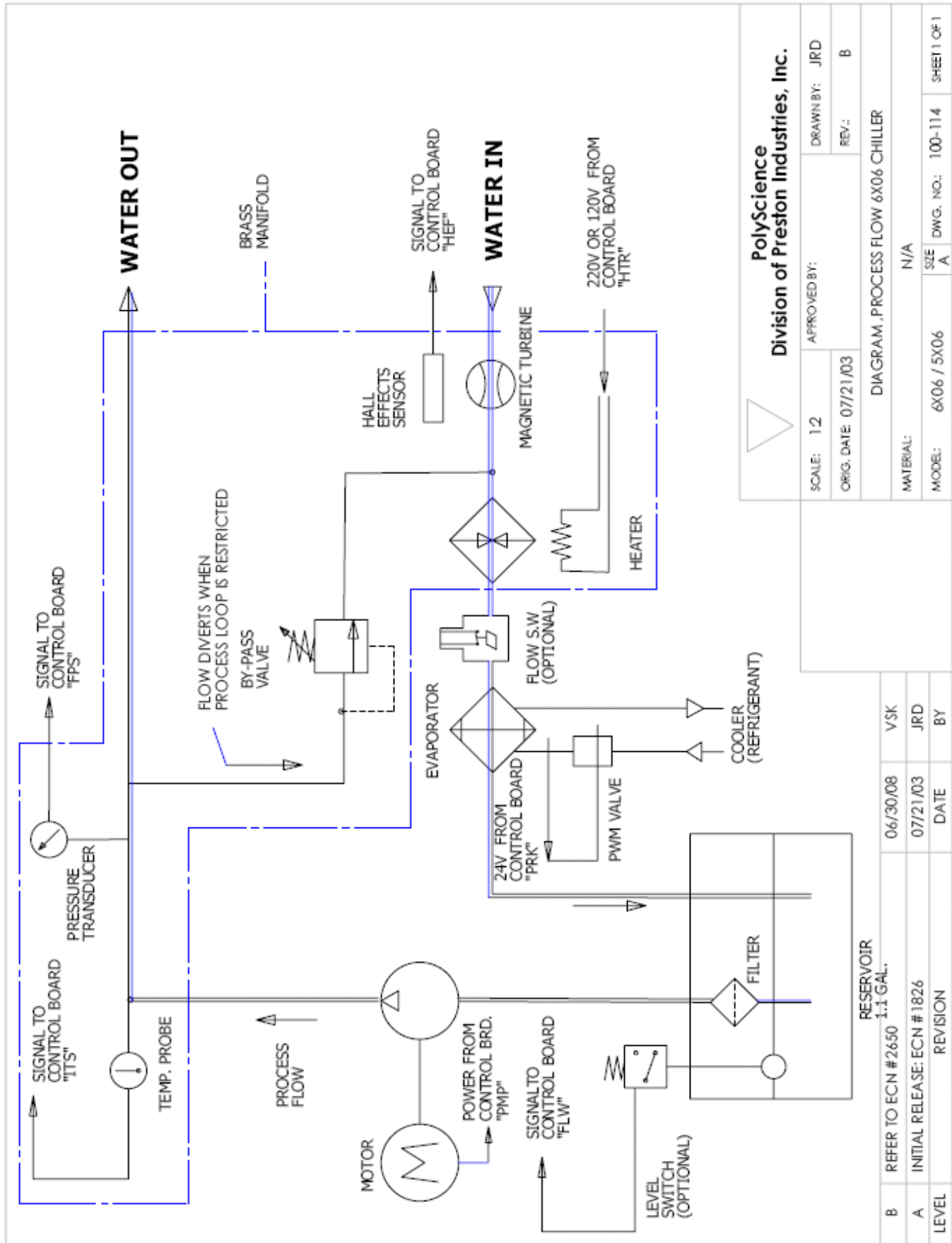

## Replacement Parts

<b>All 1/4-HP Units</b>	<b>120V, 60Hz</b>	<b>240V, 50Hz</b>
Condensing unit, 1/4 HP	750-157	750-158
Magnetic Drive Pump (models without heat)	525-551	525-552
Magnetic Drive Pump (models with heat)	525-553	525-554
Positive Displacement Motor (all models)	215-535	215-529
Positive Displacement Pump (all models)	215-105	215-105
Turbine Pump (all models)	215-656	215-305
Circuit Breaker	215-330	215-330
Fan blade	215-748	215-031
Fan motor	215-704	215-030

<b>All 1/3-HP Units</b>	<b>120V, 60Hz</b>	<b>240V, 50Hz</b>
Condensing unit, 1/3 HP	750-306	750-189
Magnetic Drive Pump (models without heat)	525-551	525-552
Magnetic Drive Pump (models with heat)	525-553	525-554
Positive Displacement Motor (all models)	215-535	215-529
Positive Displacement Pump (all models)	215-105	215-105
Turbine Pump (all models)	215-656	215-305
Circuit Breaker (magnetic drive pump models)	215-330	215-330
Circuit Breaker (positive displacement and turbine pump models)	215-388	215-330
Fan motor	215-730	215-724

<b>All 1/2-HP Units</b>	<b>120V, 60Hz</b>	<b>240V, 50Hz</b>
Condensing unit, 1/2 HP	750-308	750-309
Magnetic Drive Pump (models without heat)	525-551	525-552
Magnetic Drive Pump (models with heat)	525-553	525-554
Positive Displacement Motor (all models)	215-535	215-529
Positive Displacement Pump (all models)	215-105	215-105
Turbine Pump (all models)	215-656	215-305
Circuit Breaker	215-388	215-330
Fan blade	215-767	215-769
Fan motor	215-768	215-708

<b>All 3/4-HP Units</b>	<b>208-230V, 60Hz</b>	<b>240V, 50Hz</b>
Condensing unit, ¾ HP	750-304	750-303
Magnetic Drive Pump (models without heat)	525-552	525-552
Magnetic Drive Pump (models with heat)	525-554	525-554
Positive Displacement Motor (all models)	215-217	215-217
Positive Displacement Pump (all models)	215-106	215-106
Turbine Pump (all models)	215-305	215-305
Circuit Breaker	215-330	215-330
Fan Blade	215-631	215-631
Fan Motor	525-731	525-731
Sensor, LM35 5/16" clip	200-503	200-503
Sensor, LM35 ½" clip	200-504	200-504

<b>All 1-HP Units</b>	<b>208-230V, 60HZ</b>	<b>240V, 50HZ</b>
Compressor, 1 HP	750-304	750-303
Magnetic Drive Pump (for models without heat)	525-552	525-552
Magnetic Drive Pump (for models with heat)	525-554	525-554
Positive Displacement Motor (all models)	215-217	215-217
Positive Displacement Pump (all models)	215-106	215-106
Turbine Pump (all models)	215-305	215-305
Circuit Breaker	215-330	215-330
Fan blade	215-631	215-631
Fan motor	525-731	525-731
Sensor, LM35 5/16" clip	200-503	200-503
Sensor, LM35 ½" clip	200-504	200-504

<b>Additional Parts</b>	
Operator's Manual	110-811
Tubing Adapter Kit	510-288
Air Filter	750-855
Fluid Filter	565-102
Flow Indicator	330-571
Reservoir Cap	300-460
Reservoir Spill Cup	300-459

## RS232 Communications

**Serial Connector** — A 9-pin D-connector is provided on the back panel of the Chiller for RS232 data communication. A serial cable that uses only the following pins should be used to connect the Chiller to the computer:

- Pin #2 — data read (data from computer)
- Pin #3 — data transmit (data to computer)
- Pin #5 — signal ground

**RS232 Protocol** — The controller uses the following RS232 protocol:

- Data bits — 8
- Parity — none
- Stop bits — 1
- Flow control — none
- Baud rate — selectable (Chiller and PC baud rates must match).

**Communications Commands** — Commands must be entered in the exact format shown. Do not send a [LF] (line feed) after the [CR] (character return). Be sure to follow character case exactly. A response followed by an exclamation point (!) indicates that a command was executed correctly. A question mark (?) indicates that the Chiller could not execute the command (either because it was in an improper format or the values were outside the allowable range). A response must be received from the Chiller before another command can be sent. All responses are terminated with a single [CR].

Command Description	Command Format	Values	Return Message
Set command echo	SEi[CR]	Echo: i = 1 No Echo: i = 0	![CR]
Set on / off	SOi[CR]	On: i = 1 Off: i = 0	![CR]
Set set point	SSxxx[CR]	x = ASCII digit	![CR]
Read set point temperature	RS[CR]		+xxx.x[CR] or - xxx.x[CR]
Read temperature	RT[CR]		+xxx.x[CR] or - xxx.x[CR]
Read probe 1 temperature	R1[CR]		+xxx.x[CR] or - xxx.x[CR]
Read probe 2 temperature	R2[CR]		+xxx.x[CR] or - xxx.x[CR]
Read temperature units	RU[CR]	C or F	C[CR] or F[CR]
Read status	RW[CR]	1 = Run 0 = Standby	1[CR] or 0[CR]
Read pressure in PSI	RP[CR]		+ xxx.x[CR]
Read pressure in kPa	RK[CR]		+ xxx.x[CR]
Read flow in GPM	RG[CR]		+ xxx.x[CR]
Read flow in LPM	RL[CR]		+ xxx.x[CR]
Read remote control voltage	RC[CR]		+ xxx.x[CR]
Read line voltage	RV[CR]		+ xxx.x[CR]
Read compressor discharge temperature (°C)	RH[CR]		+xxx.x[CR] or - xxx.x[CR]
Read remote probe temperature	RR[CR]		+xxx.x[CR] or - xxx.x[CR]
Read ambient temperature on PCB	RA[CR]		+xxx.x[CR] or - xxx.x[CR]

Command Description	Command Format	Values	Return Message
Read float switch status	RX[CR]	0 = fluid level is ok 1 = fluid level is low	0[CR] or 1[CR]
Read evaporator inlet temperature (°C)	REI[CR]		+xxx.x[CR] or - xxx.x[CR]
Read evaporator outlet temperature (°C)	REO[CR]		+xxx.x[CR] or - xxx.x[CR]
Read compressor outlet temperature (°C)	RUT[CR]		+xxx.x[CR] or - xxx.x[CR]
Read fault status (see <i>Display, Alarm and Error Messages</i> )	RF[CR]	00 = System OK 18 = Standby mode 02 – 17 = Fault	![CR]
Read all PID parameters	R?[CR]		SET=xx.x, P1=xx.x, P2=xxx.x, Amb=xxx EVi=xx.x, EVo=xx.x, CHS=xxx.x HEAT P=36.0, I=60.0, D=1.0 REFRIDGE P=4.0, I=60.0, D=1.0 EVAP P=0.1, I=3.0, D=2.5 COOL P=1.0, I=20.0, D=1.0 FAN P=1.0, I=150.0, D=1.0 FAULT=00 VERSION=F5-022 ?R?

## Equipment Disposal (WEEE Directive)



or



This equipment is marked with the crossed out wheeled bin symbol to indicate it is covered by the Waste Electrical and Electronic Equipment (WEEE) Directive and is not to be disposed of as unsorted municipal waste. **Any products marked with this symbol must be collected separately, according to the regulatory guidelines in your area.**

It is your responsibility to correctly dispose of this equipment at lifecycle-end by handing it over to an authorized facility for separate collection and recycling. It is also your responsibility to decontaminate the equipment in case of biological, chemical and/or radiological contamination, so as to protect the persons involved in the disposal and recycling of the equipment from health hazards. By doing so, you will help to conserve natural and environmental resources and you will ensure that your equipment is recycled in a manner that protects human health.

Requirements for waste collection, reuse, recycling, and recovery programs vary by regulatory authority at your location. Contact your local responsible body (e.g., your laboratory manager) or authorized representative for information regarding applicable disposal regulations.

## Service and Technical Support

If you have followed the troubleshooting steps and your Recirculating Chiller fails to operate properly, contact the supplier from whom the unit was purchased. Have the following information available for the customer service person:

- Model, Serial Number, and Voltage (from back panel)
- Date of purchase and your purchase order number
- Suppliers' order number or invoice number
- A summary of your problem

## Warranty

The manufacturer agrees to correct for the original user of the product, either by repair (using new or refurbished parts), or at the manufacturer's election, by replacement (with a new or refurbished product), any defects in material or workmanship which develop during the warranty period. The standard warranty is twenty-four (24) months after delivery of the product. In the event of replacement, the replacement unit will be warranted for the remainder of the original warranty period or ninety (90) days, whichever is longer. For purposes of this limited warranty, "refurbished" means a product or part that has been returned to its original specifications. In the event of a defect, these are your exclusive remedies.

If the product should require service, contact the manufacturer's/supplier's office for instructions. When return of the product is necessary, a return authorization number is assigned and the product should be shipped, transportation charges pre-paid, in either its original packaging or packaging affording an equal degree of protection to the indicated service center. To insure prompt handling, the return authorization number must be placed on the outside of the package. A detailed explanation of the defect should be enclosed with the item.

The warranty shall not apply if the defect or malfunction was caused by accident, neglect, unreasonable use, improper service, acts of God, modification by any party other than the manufacturer, or other causes not arising out of defects in material or workmanship.

**EXCLUSION OF IMPLIED WARRANTIES.** THERE ARE NO WARRANTIES, EXPRESSED OR IMPLIED, INCLUDING, BUT NOT LIMITED TO, THOSE OF MERCHANTABILITY OR FITNESS FOR A PARTICULAR PURPOSE WHICH EXTEND BEYOND THE DESCRIPTION AND PERIOD AS STATED IN THE OPERATOR'S MANUAL INCLUDED WITH EACH PRODUCT.

**LIMITATION ON DAMAGES.** THE MANUFACTURER'S SOLE OBLIGATION UNDER THE WARRANTY IS LIMITED TO THE REPAIR OR REPLACEMENT OF A DEFECTIVE PRODUCT AND THE MANUFACTURER SHALL NOT, IN ANY EVENT, BE LIABLE FOR ANY INCIDENTAL OR CONSEQUENTIAL DAMAGES OF ANY KIND RESULTING FROM USE OR POSSESSION OF THIS PRODUCT.

Some states do not allow: (A) limitations on how long an implied warranty lasts; or (B) the exclusion or limitation of incidental or consequential damages, so the above limitations or exclusions may not apply to you. This warranty gives you specific legal rights and you may have other rights that vary from state to state.

**Air Cooled:**

6750TG5CAR93	1234950 CHILLER, 3/4HP,240/50,COHERENT
6760TG5CAR94	1234960 CHILLER, 3/4HP,208/230/60 or 200/50,COHERENT

**Water Cooled:**

6750TG5CWR93	1235050 CHILLER, 3/4HP,240/50,WC COHERENT
6760TG5CWR94	1235060 CHILLER, 3/4HP,208/230/60 or 200/50, WC COHERENT

**FILTER CHANGE REMINDER FOR AIR AND WATER FILTERS**

The filter timer counts days. For this purpose a day is defined as 24 hours. For each 24 running hours the timer is decremented by 1. For each 48 hours in standby the timer is decremented by 1. When power is off the timer holds the count that was present when power was removed.

To view the number of days remaining before change filter comes on press and hold the left hand key until HL appears scroll down. Continue to scroll down the menu until PC appears. Press one more time and CF will appear with days remaining shown on the right hand display.

To change the default value press and hold the units/menu select key until HL appears then scroll down to CF. Then press and hold while di appears. Press knob and FAC will appear. Release knob. Press left key and tF will appear Turn knob to set the desired days for the change filter reminder. Press knob when done.

When the filter days have elapsed an audible alarm will sound and the right hand display will alternately flash CHA FIL. To silence the audible press left hand key while it is sounding. CHA FIL will continue to flash. To cancel the message after checking the air and water filters go down the menu to PC then press and hold. CF will appear. While holding the key, press and release knob. Filter timer will go back to its default value.

To demonstrate the filter reminder you can force the days remaining to 0. Turn circuit breaker off then press and hold the power key and the knob while turning the circuit breaker on. Go back to previous step to cancel.

**FLUID LEVEL ALARMS**

This chiller has a dual float switch. When the fluid level falls below the first float the display will show EFL and the alarm will sound. The chiller will continue to run. If the fluid level falls below the second float the display will show EF2 and the chiller will turn off in 30 seconds if fluid is not added.

**FILLING THE CHILLER WITH FLUID**

- 1) Before the chiller is turned on, add pre-mixed Optishield to the water tank and top off on the fill port.
- 2) Run the chiller with laser head connected, run the chiller for a minute, then turn the chiller off. It may beep with low water level alarm.
- 3) Add more pre-mixed Optishield to the water tank and top off on the fill port.
- 4) Turn on the chiller again.

**CAUTION**

**Do not top off the reservoir with the chiller running. It could overflow when turned off.**

## Service Procedure for PolyScience Chillers

- 1) Check the part number and operating voltage/current requirement labeled on the chiller crate, make sure the chiller is the right model ordered, and with appropriate AC power available in the customer's lab. With a given model such as 240VAC/50Hz, **the maximum and minimum allow operation range will be +/-10% of the specified AC** (216-264VAC for the 240VAC/50Hz model for instance). Contact the factory for assistance if a wrong chiller is ordered or delivered.
- 2) **Never set the chiller temperature lower than the dew point** in the lab which may result in condensation and laser contamination/damage. Since the dew point depends on ambient temperature and relative humidity, each application must be analyzed. If the optimum chiller temperature from the laser test report is too close or already lower than the dew point in the lab, please resolve the lab condition prior to the chiller/laser operation.
- 3) Uncrate the chiller. Reference the pictures below to slide the chiller to outside. Enclosed find the user manual, chiller test report, power cords, a spare water filter and 4x 1 gallon Optishield pre-mixed (10% Optishield/90% distilled water) on the back of the crate. Double check the chiller model as labeled on the chiller back panel to ensure it matches with that on the crate. Avoid direct sun light to the Optishield for storage. Save the crate for future possible shipment of the chiller.



- 4) Please read and emphasize the following warning label to the customer.





## WARNING!

### ATTENTION COHERENT CUSTOMER

Read the following carefully otherwise any laser & chiller damage will not be covered by the warranty:

- Chiller must be operated within  $\pm 10\%$  of the nominal AC voltage at the specified line frequency.
- Optishield pre-mix 400315 is the only cooling fluid authorized by Coherent to use with this chiller.
- Any unauthorized chiller hose modification is prohibited which may lead to insufficient cooling.

- 5) Find appropriate power outlet in the customer lab to satisfy the rated AC voltage/current. 2 power cords (1 with connector, the other one no connector) are supplied with the chiller. If needed, make or buy a power cord with appropriate plug.

**Note:** The PolyScience chillers need  $\sim 220\text{VAC}/50\text{Hz}$  or  $\sim 220\text{VAC}/60\text{Hz}$  to operate. In case only  $\sim 110\text{VAC}$  utility is available, the alternative option is to prepare a step-up transformer. When connecting the transformer, ensure the AC power entry switch is correctly set (typically there is switch for  $\sim 110\text{VAC}$  or  $\sim 220\text{VAC}$  input), and ensure the AC output is at  $\sim 220\text{VAC}$  level before plug into the chiller; also make sure the  $\sim 110\text{VAC}$  power outlet and the step-up transformer are capable to yield the rated AC power for the chiller plus  $\sim 50\%$  margin for safety.

- 6) Ensure the cooling fluid not contaminated with the previous chiller/laser combination. If contamination is observed, take appropriate action to clean. Drain the cooling fluid thoroughly from the laser head to minimize the concern of the cooling fluid incompatibility.

**Warning:** In no case should vinegar, acid, Chloramine-T (corrosive to copper) or other like detergents be used!

- 7) For Evolution-HE, directly connect the chiller hoses to the laser head; **for Evolution-15/30/45, open the top cover of the chiller, back off (turn counter-clockwise) the inside flow valve by 2 turns to prevent the excessive pressure on the laser head**, then connect the chiller hoses.

**Note:** if the flow valve adjustment step is missed for Evolution-15/30/45 lasers, a fault code 08 may display and then shut off the chiller, which means the pressure exceeds the high pressure limit setting (default to 70PSI). Open the top cover and adjust the valve counter-clockwise to clear this fault. Reference the chiller Operator's manual for other fault information if in case it happens.

- 8) Open the chiller filter housing to ensure a filter cartridge sits properly; the filter cartridge should be brand new w/o loose material on either the inner or the outer wall, and no stagnant water inside the chiller housing.
- 9) For water-cooled chillers, connect the facility water hoses. **Ensure the facility water valve is open and the pressure is set appropriately (minimum 50PSI and do NOT exceed 100PSI).**

**Note:**

a) If the facility water is not connected or the facility water valve is not open before turning on the chiller, the chiller may get overheated with some burnt smells, and the water temperature may go up rather than cool to the set temperature.

b) Since the facility water loop is sealed, the pressure reading on the facility water supply or return port may not necessarily mean there is flow (for instance, even when the facility water valve is closed). When there is flow, there should be a pressure drop between the supply and the return port of the facility water loop.

- 10) Connect the water hose kit with SUPPLY port to the bottom and RETRUN from the top at the Evolution laser head. Ensure there is no open loop in the Evolution as well as the water line going to the amplifier. Then follow the fill process below:

a) Before the chiller is turned on, add pre-mixed Optishield to the water tank and top off on the fill port.

b) Run the chiller with laser head connected, run the chiller for a minute, then turn the chiller off. It may beep with low water level alarm.

c) Add more pre-mixed Optishield to the water tank and top off on the fill port.

d) Turn on the chiller again.

**Caution:** Do not top off the reservoir with the chiller running. It could overflow when the chiller is turned off.

- 11) Chiller is ready for use. Record the date (as start date of the new filter and new cooling fluid for future reference on regular preventive maintenance), the chiller pressure, the flow rate when connected with the given laser system, by toggling the menu/unit button on the left side. For reference, below is the recommendation of the chiller pressure/flow rate with different laser head. Adjust the flow valve inside if necessary. The flow display has a time lag when adjusting the bypass valve. Wait ~30 seconds for the flow reading to respond.



## Preventive Maintenance

Cooling Fluid Type \_\_\_\_\_ Date \_\_\_\_\_

By \_\_\_\_\_ Next Due NO LATER THAN \_\_\_\_\_

## Chiller Settings

Pressure \_\_\_\_\_ PSI Set Temp \_\_\_\_\_ °C/°F

\_\_\_\_\_ KPa X 100

Flow Rate \_\_\_\_\_ GPM

\_\_\_\_\_ LPM

1255902

Rev. AA

## Recommended pressure/flow rate settings for different Evolution lasers:

	For Evolution-15	For Evolution-30/45	For Evolution-HE
<b>Pressure Reading (PSI)</b>	30-40	40-50	55-65
<b>Flow Rate Reading (GPM)</b>	~1.8	~2.3	~3.8

- 12) Every time when the chiller is turned on, it shows “CF” on the left side and “xxd” on the right side, where “xx” represents how many days left before the due date of the next comprehensive preventive maintenance (including draining the cooling fluid, changing the filter, and refilling the tank etc.). The default interval is set as 90 days – note the 90 days are not the actual calendar days. For each 24 chiller running hours the timer is decremented by 1; for each 48 hours in stand-by mode the is decremented by 1. As a consequence, for OEM customers who run the chiller 7/24, the timer will beep and flash every 90 days of continuous chiller operation. For scientific customers, if assuming 8 hours a day and 5 days a week in operation, the timer will beep and flash every 145 calendar days.

The chiller will remind and beep with “CH FIL” message when the timer is counted down to 0 day. Follow the instruction below for the preventive maintenance. Also take this opportunity to check and clean the air filter (which is washable) on the front panel for air-cooled chillers.

## Reminder for all Evolution pump lasers

To prevent corrosion and algae growth in the Evolution closed-loop cooling water system, it is important to **do preventive maintenance by follow the filter change reminder, or at least every six months** (in case the chiller is turned off, the filter change reminder will not count), **whichever comes first**.

Directions:

- Drain used fluid from the chiller.
- Flush chiller with **clean steam-distilled water** for 30 minutes.
- Change particle filter.
- Clean up the condenser, the air vents and reusable air filter per Operator's Manual.
- Follow the fluid fill process
  - Note: only pre-mixed Optishield is authorized by Coherent for PolyScience chillers. You can order pre-mixed OPTISHIELD 1-gallon size from Coherent (Part #: 1166400)
  - Warning: OPTISHIELD may be harmful if swallowed, inhaled, or absorbed through the skin or eyes.
- The system is ready for use.

If the laser is used frequently, leaving the chiller on between periods of use will help prevent algae growth. If the chiller is not intended to use for more than 2 weeks, it is recommended to have it completely drained (no residual fluid in the tank, or in the filter housing).