

PRIORITY	
<input type="checkbox"/> Mandatory	<input type="checkbox"/> As Required
<input type="checkbox"/> Next Visit	<input type="checkbox"/> At Installation
<input checked="" type="checkbox"/> Information	

FIELD SERVICE BULLETIN

DATE ISSUED: 8/27/09

NUMBER: 585

Chameleon Vision Sealed Cavity

PRODUCT All Chameleon Vision systems

PURPOSE To inform the field-service network and our representatives that the Chameleon Vision is a sealed unit.

DESCRIPTION The Chameleon Vision cavity and pre-comp section is hermetically sealed in the clean-room and should never be opened in the field. Opening the cavity or pre-comp section will result in contamination and necessitate a major repair with mechanical strip-down for complete re-cleaning. The outer-lid can be removed to gain access to the Chameleon headboard and Verdi pump-board as shown in Fig1. The internal covers should never be removed as highlighted in Fig2.


Fig 1 – outer cover



Fig 2 – Internal covers

Chameleon Cavity
Sealed Area

Vision Pre-comp Cavity
Sealed Area

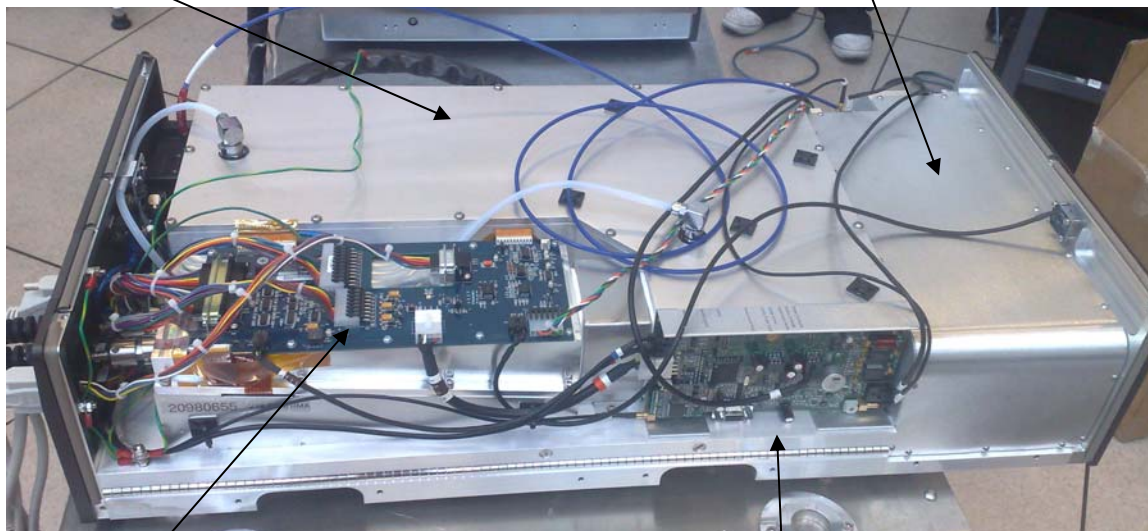

Pump Board

Head Board

ACTION

PARTS

PSE Initials IM