

## Typical Applications

- ✓ Astrophotography
- ✓ Sky Surveys
- ✓ Radiology



Apogee Aspen Series

## System Features<sup>\*1</sup>

- **High Resolution Sensor**  
9.3 Megapixel sensor with 12  $\mu\text{m}$  pixels delivers an exceptionally large field of view with high resolution.
- **Programmable TE cooling down to 60°C below ambient**  
Ideal for detection of weak chemiluminescence or astronomy images, enabling long exposure acquisitions with optimized signal to noise ratio.
- **Ethernet interface with built-in web server**  
Remote access and control over the Internet, via standard web browser.
- **USB 2.0 interface**  
Direct 'Plug and Play' simplicity of USB 2.0.
- **16-Bit digitization**  
High photometric accuracy.
- **High longevity shutter**  
Shutter during readout and take dark reference frames - 58 mm.  
Specified for >5 million cycles.
- **Programmable I/O port**  
Synchronization with intricate experimental set-ups.
- **Remote Triggering**  
LVTTTL input allows exposure to start within 25 microseconds of the rising edge of the trigger.
- **Focusing mode**  
Faster readout option, ideal for focus optimisation.
- **Precision locking filter wheels optional**  
Choose from a range of Apogee family filter wheels with up to 17 positions.

## Apogee Aspen CG9000: Compact, 9.3 Megapixel CCD

Ideally suited to challenging astronomy and physical science imaging applications, the Apogee Aspen family offers a range of popular full frame and interline CCD sensors, within a camera platform that is designed to push performance. Deep thermoelectric cooling ensures optimal sensitivity for long exposure applications. The simple convenience of a USB 2.0 interface is accompanied by an Ethernet network interface with a built-in web server. The Apogee Aspen also utilizes a new extremely high reliability shutter, specified for > 5 million shutter cycles.

The Aspen CG9000 uses a very large format 9.3 megapixel full frame sensor with 12  $\mu\text{m}$  pixels and anti-blooming gates. Cooling down to 60°C below ambient results in a low dark current contribution. These features make the Aspen CG9000 an ideal solution for applications requiring both a large field of view and optimal signal to noise ratio, such as astrophotography, sky surveys and radiology.

## Specifications Summary<sup>\*1</sup>

<b>Array Size (pixels)</b>	3056 x 3056 (9.3 Megapixel)
<b>Pixel Size</b>	12 x 12 $\mu\text{m}$
<b>Sensor Size</b>	36.7 x 36.7 mm (1345 mm <sup>2</sup> ) 51.6 mm diagonal
<b>Pixel Well Depth (typical)</b>	90,000 e <sup>-</sup>
<b>Dark Current<sup>*2</sup></b>	0.0116 e <sup>-</sup> /pixel/sec
<b>Read Noise<sup>*3</sup></b>	9.4 e <sup>-</sup> (RMS @ 0.87 Mhz)
<b>Maximum Dynamic Range</b>	79.6 dB (9574:1)
<b>Quantum Efficiency</b>	64% @550 nm 37% @400 nm

## SPECIFICATIONS


### Technical Specifications<sup>\*\*</sup>

Sensor Type	KAF-09000 (ON Semiconductor)
Active pixels	3056 x 3056 W x H (9.3 Megapixel)
Sensor Size	36.7 x 36.7 mm (1345 mm <sup>2</sup> ) 51.9 mm diagonal
Pixel Size	12 x 12 $\mu$ m
Pixel Well Depth	90,000 e <sup>-</sup>
Read Noise <sup>**3</sup>	9.4 e <sup>-</sup> (RMS @ 0.87 MHz)
Pixel Binning	1 x 1 to 8 x 3056 on chip
Quantum Efficiency <sup>**4</sup>	64% @550 nm 37 % @400 nm
Cooling	Maximum cooling up to 60°C below ambient temperature; -35°C at 25°C ambient Thermoelectric cooler with forced air.
Temperature Stability	+/- 0.1°C
Dark Current <sup>**3</sup>	0.0116 e <sup>-</sup> /pixel/sec
Blemish Specification	Grade S as per sensor manufacturer definition
Anti-blooming factor	>100x
Maximum Dynamic Range	79.6 dB (9574:1)
Linearity	Better than 99%
Frame Rate (fps) <sup>**5</sup>	0.09 Full frame (@0.87 MHz) 0.43 Full frame (@4.17 MHz, focusing mode)
Frame Sizes	Full frame, sub-frame
Digital Resolution	16-bit
Camera Window	UV-grade fused silica

### General Specifications

Interface Options	USB 2.0 Ethernet: Network interface with built-in web server, up to 2 MHz throughput
Remote Triggering	LVTTTL trigger input, expose strobe output
Peripheral communications	8 pin mini-DIN I/O connector
Image Sequencing	1 to 65535 image sequences under software control
Exposure Time	100 milliseconds to 183 minutes (2.56 microsecond increments)
Operating System Support	Windows, Linux

## Quantum Efficiency (QE) Curve\*\*


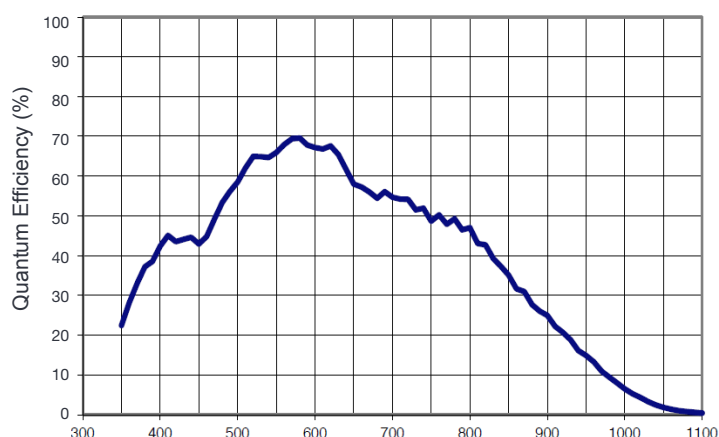


## Size of CCD Imaging Area

36.7 x 36.7 mm



## Mechanical Drawings

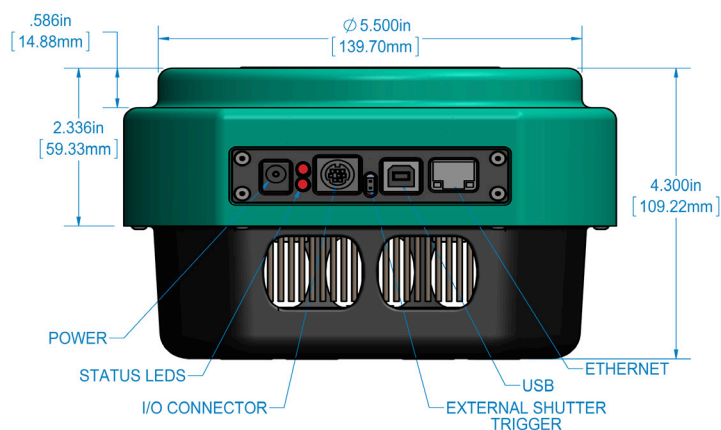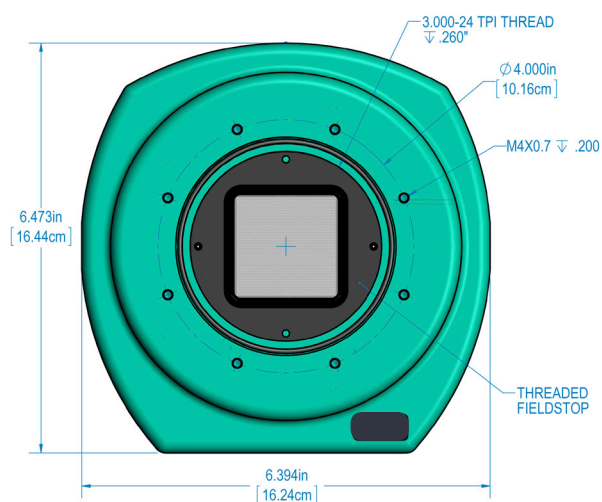


## Mechanical Specifications

Camera Housing	Aluminum, hard anodized (G07)
Camera Head Size	6.5" x 6.4" x 4.3" [16.4 x 16.2 x 10.9 cm]
Back Focal Distance	1.013" (2.57 cm) [optical]
Mounting	2.25" Aperture, 3" Thread Option and Threaded Fieldstop
Shutter	58 mm shutter (specified for >5 million cycles)
Weight	3.1 lb. (1.4 kg)

## CREATING THE OPTIMUM PRODUCT FOR YOU

How to customize the Apogee Aspen CG9000:

### Step 1: Select your camera type

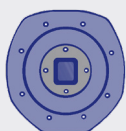


Camera

Description	Part Code
Apogee Aspen CG9000 9.3 Megapixel CCD camera with grade S sensor and 58 mm Shutter	CG9000-S-G07-S58



### Step 2: Please indicate which adapters and accessories are required

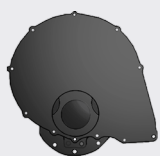

Adapters &  
Accessories

A wide range of mounting adapters and accessory options are available for the Aspen. Please refer to the links below on the Andor website for further information on filter wheels, filters and adapters.

#### Filter Wheels

Filter wheels available with up to 17 filter positions.

Please refer to [Apogee Filter Wheels](#)

#### Filters

A comprehensive selection of Astrodon filters and filters are available to complement your selected filter wheel

Please refer to [Apogee Filters](#)

#### Lens Adapters and flanges

Select the required camera mounting option for your application, from our range of lens, telescope and slip-fit faceplate adapters.

Please refer to [Apogee Adapter Matrix](#)



### Step 3: Please indicate which software you require



Software

The Apogee Aspen also requires at least one of the following software options:

Description	Ordering Information
Windows SDK for Apogee	Please download from the <a href="#">Apogee Downloads Page</a>
ASCOM Camera and Filter Wheel Driver	Please download from the <a href="#">Apogee Downloads Page</a>
Linux Driver CD	40053
Maxim DL Pro Software CD	40054
MicroManager	Please see <a href="https://micro-manager.org/wiki/Apogee">https://micro-manager.org/wiki/Apogee</a>

# Order Today

Need more information? At Andor we are committed to finding the correct solution for you. With a dedicated team of technical advisors, we are able to offer you one-to-one guidance and technical support on all Andor products. For a full listing of our local sales offices, please see: [andor.com/contact](http://andor.com/contact)



Our regional headquarters are:

## Europe

Belfast, Northern Ireland  
Phone +44 (28) 9023 7126  
Fax +44 (28) 9031 0792

## Japan

Tokyo  
Phone +81 (3) 6732 8968  
Fax +81 (3) 6732 8939

## North America

Connecticut, USA  
Phone +1 (860) 290 9211  
Fax +1 (860) 290 9566

## China

Beijing  
Phone +86 (10) 8271 9066  
Fax +86 (10) 8271 9055



**Front page image:** Melotte 15 in the Heart Nebula, courtesy of Dan Goldman. Check out other astounding images with Apogee cameras on the Andor website.

## Footnotes

1. Figures are typical unless stated otherwise
2. At minimum temperature
3. Readout noise is for the entire system. It is a combination of sensor readout noise and A/D noise.
4. Quantum efficiency of the sensor at 25°C, as supplied by the sensor manufacturer.
5. Assumes internal trigger mode of operation and minimum exposure time.



### PC Requirements

- 3.0 GHz single core or 2.4 GHz multi core processor
- 2 GB RAM
- 100 MB free hard disc to install software (at least 1GB recommended for data spooling)
- USB 2.0 High Speed Host Controller capable of a sustained rate of 40MB/s
- Windows (XP, Vista, 7 and 8) or Linux

### Operating and Storage Conditions

- Operating Temperature: 0 to 40°C
- Relative Humidity: < 70% (non-condensing)
- Storage Temperature: -25°C to 50°C
- Altitude up to 2000 m

### Power Requirements

- 100-240V, AC 50-60Hz, or alternate 12V input from user's source.
- 75W maximum power consumption (shutter open and cooling maximum)

