


Q-smart HE

Compact High-Energy pulsed Nd:YAG lasers
with excellent beam quality and versatility

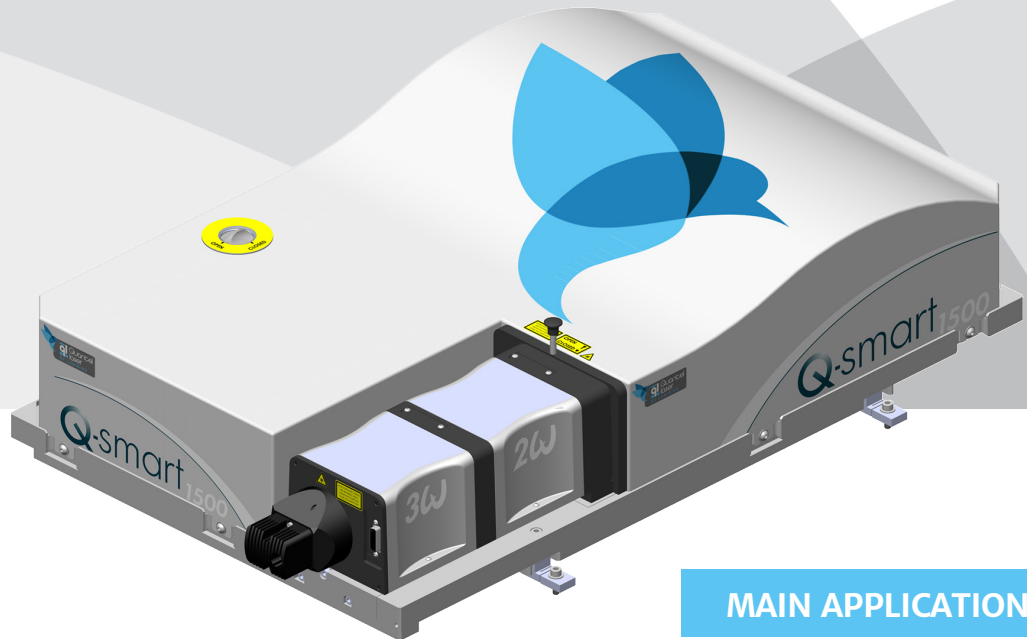


MAIN FEATURES

- Up to 2.3 J @ 1064 nm
- Robust and field proven technology
- Built to last thanks to ceramic reflectors and long flashlamp life time warranty
- Plug & play harmonic modules with automatic phase-matching
- Cables and cooling lines fully disconnectable
- Easy to use and maintain
- Optional chiller
- Intuitive GUI interface
- SLM option (Single Longitudinal Mode)


MAIN APPLICATIONS

- LiDAR
- INSTRUMENTATION
- PLD
- DYE, OPO & Ti:Sa PUMPING
- SPECTROSCOPY
- LIF
- COMBUSTION


Typical beam profiles

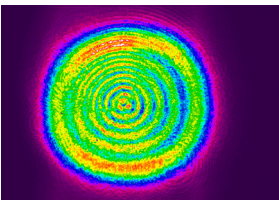

Near field 1.5 J @ 1064 nm, 10 Hz

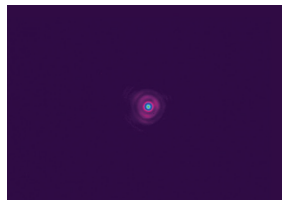

Far field 1.5 J @ 1064 nm, 10 Hz

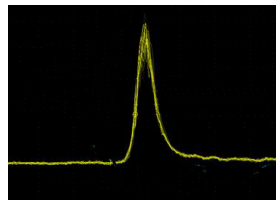

6 ns typical temporal profile
@ 1064 nm
(1 GHz oscilloscope)

www.quantel-laser.com

Many options and configurations are available.
Please contact Lumibird to find the best match for
your needs and compatibility between options.

Q-smart HE

Compact High-Energy pulsed Nd:YAG lasers with excellent beam quality and versatility




SPECIFICATIONS

		Q-smart 1200		Q-smart 1500		Q-smart 2300
Repetition rate (Hz)		10	20	10	20	10
Energy per pulse (mJ)	1064 nm	1200	1000	1500	1400	2300
	532 nm	575 / 650 ⁽¹⁾	480 / 540 ⁽¹⁾	750 / 850 ⁽¹⁾	660 / 750 ⁽¹⁾	1150/1300 ⁽¹⁾
	355 nm	280 / 350 ⁽²⁾	220 / 270 ⁽²⁾	400 / 520 ⁽²⁾	360 / 460 ⁽²⁾	620 / 850 ⁽²⁾
	266 nm	110	80	130	100	200
Pulse duration (ns) ⁽³⁾	1064 nm	5 - 10				
Beam diameter (mm) ⁽⁴⁾	1064 nm	≤ 10				≤ 12
Beam divergence (mrad) ⁽⁵⁾	1064 nm	≤ 0.5				
M ² ⁽⁶⁾	1064 nm	≤ 2				
Spatial profile @ 1064 nm ⁽⁷⁾ (fit to Gaussian)	Near field ⁽⁸⁾	≥ 0.7				
	Far field ⁽⁹⁾	≥ 0.9				
Polarization ratio (%) ⁽¹⁰⁾	1064 nm	≥ 90				≥ 80

- (1) 532 nm high energy version
 (2) 355 nm high energy version
 (3) Measured at FWHM with fast photodiode and 1 GHz oscilloscope
 (4) At the output of the laser
 (5) Full angle, at 1/e² of the peak
 (6) At 1/e² of the peak, measured by Spricon LBA FWB
 (7) Least square fit to Gaussian (perfect fit = 1)
 (8) Measured at 1 m from laser output
 (9) Measured at focal plane of a 2 m focus lens
 (10) Polarization is horizontal @ 1064, 355 & 266 nm and vertical @ 532 nm

Pulse to pulse energy stability (%) ⁽¹¹⁾	1064 nm	± 2 (0.6)
	532 nm	± 4 (1.3)
	355 nm	± 6 (2)
	266 nm	± 8 (2.6)
Power drift (%) ⁽¹²⁾	1064 nm	± 3
	532 nm	± 5
	355 nm	± 5
	266 nm	± 10
Pointing stability (μrad) ⁽¹³⁾	1064 nm	< 40
Linewidth @1064 nm (cm ⁻¹)	Standard ⁽¹⁴⁾	≤ 0.7
	SLM ⁽¹⁵⁾ option	≤ 0.005

- (11) Peak-to-peak (RMS), 100% of shots
 (12) Over 8 hours for ΔT° ≤ ± 3°C
 (13) Measured with Spiricon LBA-100, rms, on 200 pulses at the focal plane of a 2 m focus lens
 (14) Measured at FWHM with a grating spectrometer with 0.045 cm⁻¹ resolution
 (15) Measured at FWHM with a slow scan Fabry-Perot etalon, ≤ 20% energy reduction @ 1064 nm

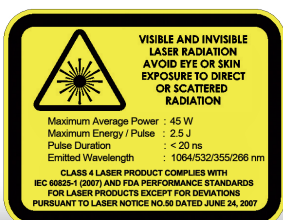


OTHER INFORMATION


Power requirements		240 V (110 V on request) 16 to 20 A, 50/60 Hz
Cooling		Water to water Optional chiller (stand alone or 19» rack)
Operating temperature		+ 18 °C to + 28 °C
Storage temperature ⁽¹⁶⁾		- 10 °C to + 50 °C
Cable length (m)		3 ⁽¹⁷⁾
Flashlamps warranty		50 million shots ⁽¹⁸⁾
Weight (kg)	Laser head	45
	Harmonic modules	2.1
	Integrated cooling & electronics	50 / 70 ⁽¹⁹⁾

- (16) System rinsed and drained with ethylene glycol/water mixture
 (17) Other lengths on request. Some losses are to be expected.
 (18) 80% of energy, or 1 year, whichever comes first
 (19) Q-smart 1200-1500 / Q-smart 2300

Laser head

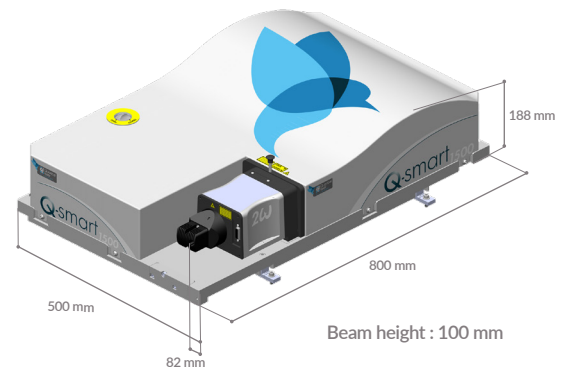

Integrated cooling & electronics

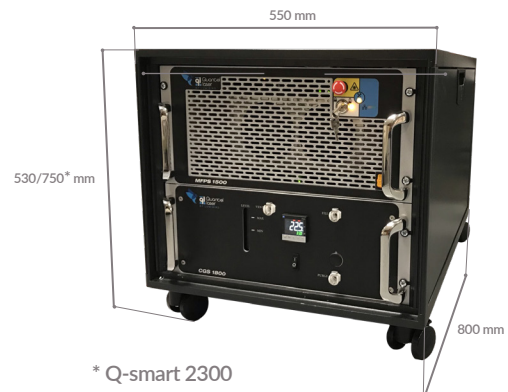

* Q-smart 2300

www.quantel-laser.com

Many options and configurations are available. Please contact Lumibird to find the best match for your needs and compatibility between options.



Lumibird has locations across the globe that are available to provide support for any product, service or inquiry. Visit www.lumibird.com to connect with any of our global sites.

